

ENOCHS Regency 60

OB/GYN EXAMINATION TABLE

Vacuum-formed Top is an 100% seamless design free of decorative embossing, cross seams and fabric folds that make cleaning difficult and inefficient. Top is padded with 1½" of high-density foam with an extra ½" on along the edges for patient comfort. Top is covered in durable Naugahyde® manufactured to be wear, stain, mold, mildew and oil resistant

Attainable Positions: OB/GYN, Lithotomy, Urology and Flat

Dimensions: Length - 52" Extends to 71"
Width - 27"
Height - 34½"

Accessories: Knee crutches, stainless steel basin, padded armboard, restraint straps, IV Pole, urology basin with screen and drain, headend accessory rails, footend accessory rails and caster base.

Paper Roll Holder is concealed at head end of table and accommodates 18" or 21" paper.

Perineal Recess is 100% seamless for infection control maintenance.

Back smoothly adjusts manually from horizontal to 87° chair position.

Pelvic Tilt - 3-positions 0°, 7°, 11° with release handles on either side of table.

Reset Button for circuit breaker system, restores power to table's outlets and warming drawer.

Positive Lock Multi-position Stirrups with length and lateral adjustment, conceal when not in use.

Warming Drawer Indicator On/Off Switch illuminates when warmer is on.

Padded **Pullout Leg Extension** works independent of irrigation basin to allow use of basin when unit is extended.

UL listed **Dual Electrical Outlets** located for quick convenient access.

4" deep **Irrigation Basin** operates independent of the leg extension for quick, easy transition.

Four 4" deep **Storage Drawers** all have rounded interior corners and edges for easy cleaning

The 4" deep **Warming Drawer** has rounded internal corners and edges for quick easy cleaning.

Table Rating: 120 V AC 60 Hz 7 A
Patient Load: 400# evenly distributed

Pullout Locking Step is covered in ribbed rubber tread for easy cleaning.

Table is UL Listed

ENOCHS is an ISO 9001:2000 and WOB Certified Company






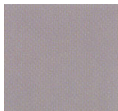
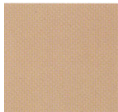









ENOCHS Vacuum-formed Top has an all steel frame and steel back design for strength, durability and long service life. The sleek efficient top features a seamless surface, moulded perineal recess and pelvic tilt positions of 0°, 7° & 11°. The pullout padded leg extension is an exclusive design that allows extension to move independently of the 4" deep irrigation basin. The separate operation allows for patient to go from extension to stirrups quickly and easily while still allowing unrestricted access to basin.

The Regency 60 has a total of 3.3 cubic feet of storage space distributed between three side and two footend drawers. ENOCHS drawers have a steel reinforced side guide system for smooth non-binding movement. Automatic stop built into the drawer eliminates accidental spillage of drawer contents.

The Regency 60's pull out step at the footend of the table is covered with easy to clean ribbed tread and locks into position.

The ELPK 1000 is located on front side of table for easy access. The assembly houses two UL LISTED Hospital Grade electrical outlets, ON/OFF switch for the warming drawer indicator light which illuminates when warming drawer is ON, and reset button which restores power to table's outlets and warming drawer if the circuit breaker system is tripped.

For additional information on the ENOCHS Regency 60 and other fine examination room furniture call Customer Service at 800-428-2305, FAX to 317-580-2944, e-mail to enochs@enochsmed.com or write to P.O. Box 50559 Indianapolis, IN 46250

Feature	ENOCHS REGENCY 60 SPECIFICATION		
<p>Color Selections</p> <p>Colors and textures represented here, are an approximate of true color under the limits of electronic viewing.</p> <p>Contact ENOCHS for actual color samples.</p> <p>Colors subject to change without notification</p>	<p>Base Colors</p>  <p>Pearl Grey</p>  <p>Cream White</p>  <p>Mocha Brown</p>	<p>Drawer Front Colors</p>  <p>Grey</p>  <p>Beige</p>  <p>Natural Oak</p>	<p>Upholstery Top Colors</p>  <p>Cambridge Blue</p>  <p>Smoke</p>  <p>English Brown</p>  <p>Windsor Wine</p>  <p>Mountain Violet</p>  <p>Desert Sand</p>
<p>Ordering Configuration</p>	 <p>If room requires drawers to open on <u>PATIENT'S</u> right. ORDER a Right Hand Table (STANDARD)</p> <p style="text-align: center;">OR</p> <p>If room requires drawers to open on <u>PATIENT'S</u> left. ORDER a LEFT Hand Table</p> 		
<p>Base Construction</p>	<ul style="list-style-type: none"> • 18 gauge steel reinforced with heavier gauge steel for strength • Fully welded base and top construction for dimensional stability • Powder coated metal surfaces for durability and long service life • Column style base with smooth surfaces for easy efficient cleaning 		
<p>Top Construction</p>	<ul style="list-style-type: none"> • Tubular steel frame with steel backing • 100% seamless vacuum-formed top design is free of cross-seams, embossing and fabric folds for quick and effective cleaning • Padded with 1½" high density foam with extra ½" contouring along edges for patient comfort • Covered in durable Naugahyde® manufactured to be wear, stain, mold, mildew and oil resistant • Upholstered sections securely attached to top frame, eliminating areas of contamination associated with removable top pads. • Top engineered for easy replacement if required 		
<p>Perineal Recess</p>	<ul style="list-style-type: none"> • Moulded recess in seat section allows greater access to patient's perineal region. during examination or procedure • Recess measures 10" wide at opening and 2" deep at center • Allows greater feeling of stable positioning for patient 		
<p>Leg extension and irrigation basin</p>	<ul style="list-style-type: none"> • Pull out padded leg extension extends table to 71" " Padboard – 19" W x 19" D • 4" deep stainless steel basin moves independent of padded leg extension for quick and easy transition • Patients with limited hip rotation may be positioned with feet on extension instead of in stirrups while physician retains full access to irrigation basin 		
<p>Manual Adjustments</p>	<ul style="list-style-type: none"> • Back adjusts to an infinite number of degrees from horizontal to 87° • 3-position Pelvic Tilt : 0°, 7°, & 11° • Leg extension adds 19" to length of table 		
<p>Stirrups</p>	<ul style="list-style-type: none"> • Multi-position with three lateral positions with quick and easy length adjustment 		
<p>Storage</p>	<p>Three 4" deep side drawers and two 4" deep foot end drawers offering 3.3 cubic feet of total storage</p> <ul style="list-style-type: none"> • High impact, rounded interior corners and edges for easy cleaning • Steel reinforced side guide mounting system for smooth non-binding movement • One-piece drawer front is covered front, back and edges with durable laminate 		
<p>Electrical</p>	<ul style="list-style-type: none"> • Dual electrical outlets located on front side of table for convenient access • Warming drawer On/Off indicator switch located on side of table for easy identification of warmer being active • Exclusive feature: circuit breaker system, automatically shuts down tables outlets if electricity is interrupted • Reset button quickly restores power to table's outlets and warming drawer 		
<p>Warmer Drawer</p>	<p>Standard feature for the Regency 60</p>		

Feature	ENOCHS REGENCY 60 SPECIFICATION
Accessories	<ul style="list-style-type: none"> • universal clamp (part# 6109) • tri-clamp (part# 6110) • headend accessory rails (part# 111-02) • footend accessory rails (part# 111-03) • knee crutches† (part# 111-11) • padded armboard†† (part# 111-05) • restraint straps (part# 4126) • IV Pole*†† (part# 6130) • Stainless Steel basin (part# 111-13) • urology pan with drain, screen & hose (part# 111-10) <p style="font-size: small; margin-top: 5px;">* requires #6109 or #6110 clamp † requires footend accessory rails #111-03 †† requires headend accessory rails #111-02</p>
Specifications	<ul style="list-style-type: none"> • Length - 52" extends to 71" • Width - 27" • Height - 34½" • 4" deep irrigation pan • Table Rating: 120 VAC 60 Hz 7A • Entire Product UL Listed
Warranty	4 – YEAR WARRANTY
Replacement Parts After Warranty Ends	ENOCHS strives to use standard commercial available parts in ENOCHS products. After warranty period ends, most parts can be purchased from ENOCHS or from local supply house which allows facility to secure the most competitive price for the parts.

ORDERING INFORMATION

ENOCHS can only be ordered through Authorized ENOCHS distributors.
 Call 800-428-2305 or e-mail enochs@enochsmed.com for the distributor in your area.

***QUICK-SHIP** Regency 60 OB/GYN Table - Right Handed** **Product Code**

Grey Drawer Fronts, Pearl Grey Base, Smoke Vacuum-formed Top	60-GRPGSM
Grey Drawer Fronts, Pearl Grey Base, Windsor Wine Vacuum-formed Top	60-GRPGWW
Grey Drawer Fronts, Pearl Grey Base, Cambridge Blue Vacuum-formed Top	60-GRPGCB

Other Popular Regency 60 Color Combinations **

Grey Drawer Fronts, Pearl Grey Base, Mountain Violet Vacuum-formed Top	60-GRPGMV
Beige Drawer Fronts, Cream White Base, English Brown Vacuum-formed Top	60-BECWEB
Beige Drawer Fronts, Cream White Base, Desert Sand Vacuum-formed Top	60-BECWDS
Beige Drawer Fronts, Cream White Base, Mountain Violet Vacuum-formed Top	60-BECWMV
Beige Drawer Fronts, Cream White Base, Windsor Wine Vacuum-formed Top	60-BECWWW
Beige Drawer Fronts, Cream White Base, Cambridge Blue Vacuum-formed Top	60-BECWCB
Grey Drawer Fronts, Mocha Brown Base, Smoke Vacuum-formed Top	60-GRMBSM

* Quick-ship products are those commonly kept in the ENOCHS inventory. Stock levels fluctuate based on orders.
 Call 800-428-2305 for availability.

** Non-inventoried product. When ordered, product will be scheduled into normal monthly production.