ENOCHS EXCEL 355

General Examination Table

PATIENT LOAD: 600 LBS.

evenly distributed on the table's main couch

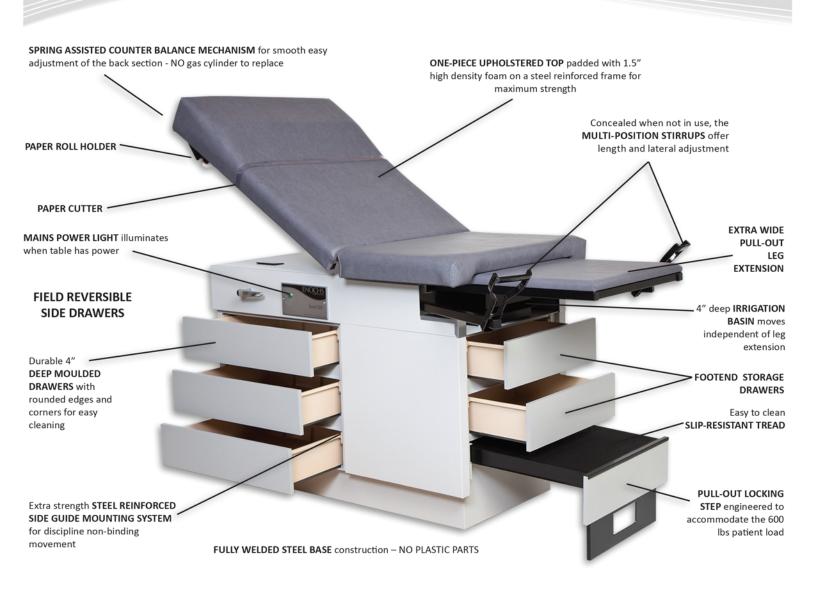


TABLE FEATURES:

- 4-YEAR Warranty
- All Steel Base Construction NO PLASTIC
- Dual Electrical Outlets, Hospital Grade, 9 amps total load
- Easy to clean One-Piece Upholstered Top, designed without cross-seams, decorative stitching or embossing
- Top constructed with a Steel Reinforced Frame for maximum strength
- Irrigation Basin moves independent of Leg Extension for fast transition
- Multi-position Stirrups with easy length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers



The extra long pull-out foot step features a steel stabilizing leg designed to accommodate a 600 pound patient load capacity. The step easily extends out from the table and locks into position providing the patient with a stable surface when getting on and off the table. The step's all steel construction offers extraordinary strength and durability. The 17.5" long surface is covered in easy to clean slip-resistant tread.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron® on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



The table's electrical receptacles feature a new tool required access door designed to restrict unauthorized use of the outlets. Door easily unlocks so accessories can be plugged into the UL Listed, hospital grade receptacles. Once in place, the door closes over the electrical cords, securing plugs in place with the locking mechanism. Receptacles have a 9 amps total maximum load.