


# ENOCHS EXCEL 356

## General Examination Table

**PATIENT LOAD: 600 LBS.**  
evenly distributed on the table's main couch

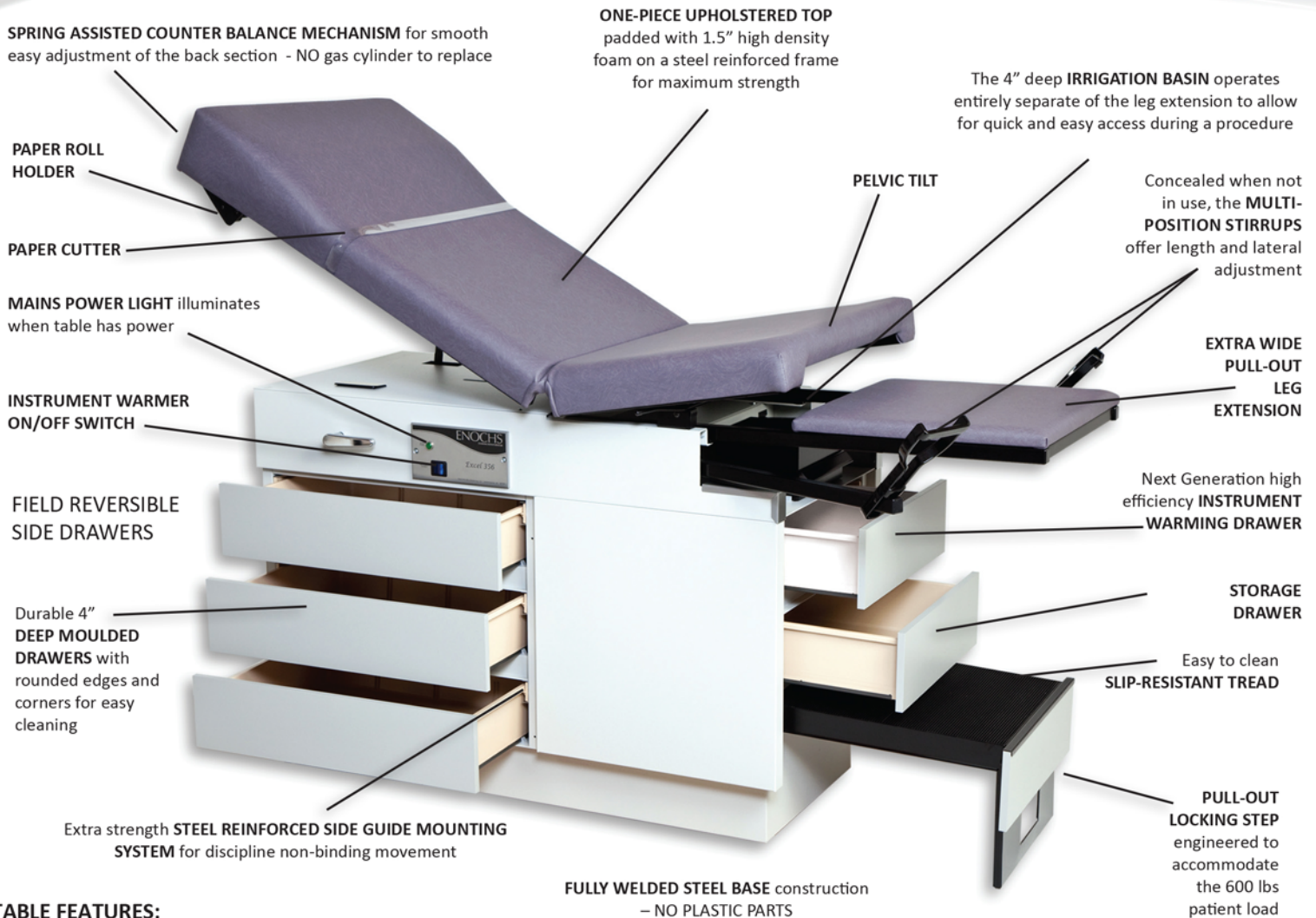

### TABLE FEATURES:

- 4-YEAR Warranty
- All welded Steel Base treated for resistance to rust
- Dual Electrical Outlets, Hospital Grade, 9 amps total load
- Instrument Warming Drawer efficiently warms drawer contents to approximate body temperature
- ON/OFF Warmer Switch with integrated Circuit Breaker - NO FUSES
- Easy to clean One-Piece Upholstered Top designed without cross-seams, decorative stitching or embossing
- 11° Pelvic Tilt with release levers on both sides of table
- Easy to adjust Stirrups with length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers



The extra long pull-out foot step features a steel stabilizing leg designed to accommodate a 600 pound patient load capacity. The step easily extends out from the table and locks into position providing the patient with a stable surface when getting on and off the table. The step's all steel construction offers extraordinary strength and durability. The 17.5" long surface is covered in easy to clean slip-resistant tread.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron® on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



The spring assisted counter balance lift mechanism (NO gas cylinder) allows for quick and convenient adjustment of the couch back to an infinite number of positions from horizontal to a full 80° chair position. The table's couch top is designed for strength and durability with a steel reinforced frame and covered in durable Naughahyde® material that provides a long, trouble-free service life.