

ENOCHS REGENCY 450

OB/GYN Examination Table

PATIENT LOAD: 600 LBS.
evenly distributed on the table's main couch

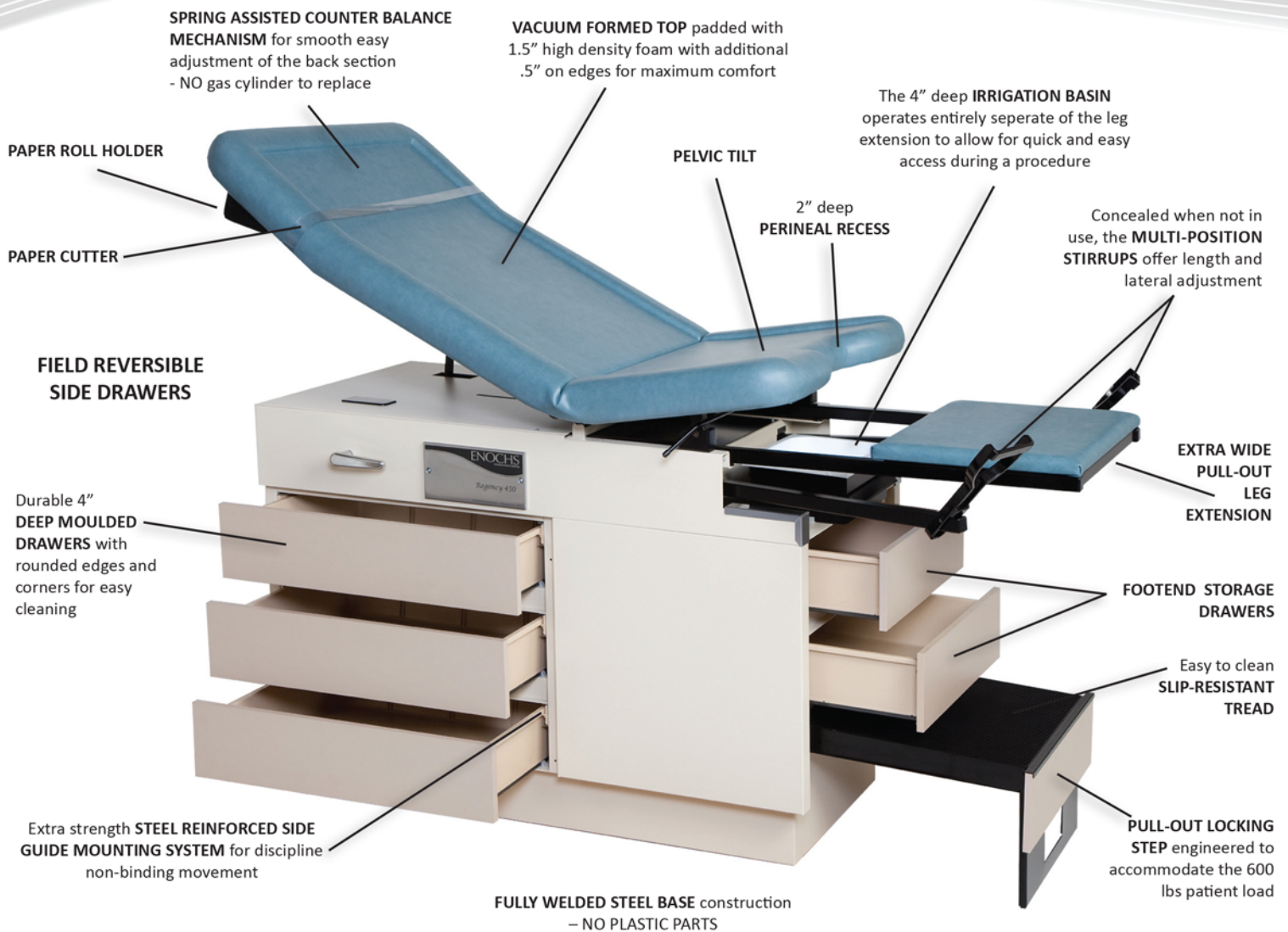


TABLE FEATURES:

- 2-YEAR Warranty
- Fully Welded all Steel Base and Top Construction - NO PLASTIC
- Spring Assisted counter balance allows smooth adjustment of couch back
- Easy to clean 100% Seamless Vacuum Formed Top
- Pelvic Tilt with a 7° and 11° adjustment
- 4" deep Irrigation Basin engineered to move independent of Leg Extension
- Easy to adjust Stirrups with length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers



The working end of the table features the independent operation between the stirrups, leg extension and irrigation basin. This efficient design gives practitioner fast and easy transition from one to another. The padded leg extension pulls smoothly out from the table to support patient's legs and can be pushed under the seat section, allowing full access to the irrigation basin and perineal recess.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron® on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



The spring assisted counter balance lift mechanism (NO gas cylinder) allows for quick and convenient adjustment of the couch back to an infinite number of positions from horizontal to a full 80° chair position. The Vacuum formed top, designed for strength and durability has an all steel back and frame (NO wooden or plastic parts). The durable Naughtaform® covered top is easy to maintain and provides a long, trouble-free service life.