ENOCHS

PRODUCT CATALOG

2021

Table of Contents

3 ENOCHS PROCEDURAL TABLES

- 4 RESOLUTE Comparison sheet
- 5 RESOLUTE 6250 Quick reference brochure
- 6 RESOLUTE 6250 Printable brochure
- 7 RESOLUTE 6250 Specifications
- 8 RESOLUTE 6220 Quick reference brochure
- 9 RESOLUTE 6220 Printable brochure
- 10 RESOLUTE 6220 Specifications
- 11 RESOLUTE 6210 Quick reference brochure
- 12 RESOLUTE 6210 Printable brochure
- 13 RESOLUTE 6210 Specifications
- 14 RESOLUTE Accessories
- 15 RESOLUTE Color selections
- 16 RESOLUTE Build order number form

ENOCHS HIGH-LOW TABLES

17 ENCORE LOW-ACCESS HIGH-LOW TABLES

- 18 ENCORE LOW-ACCESS Comparison sheet
- 19 ENCORE LOW-ACCESS 4550 Quick reference brochure
- 20 ENCORE LOW-ACCESS 4550 Printable brochure
- 21 ENCORE LOW-ACCESS 4550 Specifications
- 22 ENCORE LOW-ACCESS 4550 Build order number form
- 23 ENCORE LOW-ACCESS 4520 Quick reference brochure
- 24 ENCORE LOW-ACCESS 4520 Printable brochure
- 25 ENCORE LOW-ACCESS 4520 Specifications
- 26 ENCORE LOW-ACCESS 4520 Build order number form
- 27 ENCORE LOW-ACCESS 4510 Quick reference brochure
- 28 ENCORE LOW-ACCESS 4510 Printable brochure
- 29 ENCORE LOW-ACCESS 4510 Specifications
- 30 ENCORE LOW-ACCESS 4510 Build order number form
- 31 ENCORE LOW-ACCESS 4500 Series accessories

32 ENCORE EASY-ACCESS HIGH-LOW TABLES

- 33 ENCORE EASY-ACCESS Comparison sheet
- 34 ENCORE EASY-ACCESS 4250 Quick reference brochure
- 35 ENCORE EASY-ACCESS 4250 Printable brochure
- 36 ENCORE EASY-ACCESS 4250 Specifications
- 37 ENCORE EASY-ACCESS 4200 Quick reference brochure
- 38 ENCORE EASY-ACCESS 4200 Printable brochure
- 39 ENCORE EASY-ACCESS 4200 Specifications
- 40 ENCORE EASY-ACCESS 4200 Build order number form
- 41 ENCORE EASY-ACCESS 4200 Series accessories
- 42 ENCORE Color selections

43 ENOCHS TREATMENT TABLES

- 44 HORIZON Comparison sheet
- 45 HORIZON 2250 Quick reference brochure
- 46 HORIZON 2250 Printable brochure
- 47 HORIZON 2250 Specifications
- 48 HORIZON 2250 Color selections
- 49 HORIZON 49, 49D Color selections
- 50 HORIZON 48 Color selections
- 51 HORIZON Build order number form

52 ENOCHS FIXED-HEIGHT TABLES

- 53 REGENCY Comparison sheet
- 54 REGENCY 456 Quick reference brochure
- 55 REGENCY 456 Printable brochure
- 56 REGENCY 456 Specifications
- 57 REGENCY 450 Quick reference brochure
- 58 REGENCY 450 Printable brochure
- 59 REGENCY 450 Specifications
- 60 REGENCY Accessories
- 61 REGENCY Color selections
- 62 REGENCY Build order number form

63 ENOCHS FIXED-HEIGHT TABLES

- 64 EXCEL Comparison sheet
- 65 EXCEL 356 Quick reference brochure
- 66 EXCEL 356 Printable brochure
- 67 EXCEL 356 Specifications
- 68 EXCEL 355 Quick reference brochure
- 69 EXCEL 355 Printable brochure
- 70 EXCEL 355 Specifications
- 71 EXCEL 350 Quick reference brochure
- 72 EXCEL 350 Printable brochure
- 73 EXCEL 350 Specifications
- 74 EXCEL Accessories
- 75 EXCEL Color selections
- 76 EXCEL Build order number form

77 ENOCHS ADDITIONAL PRODUCTS

- 78 Physician Seating specifications
- 79 Physician Seating color selections
- 80 Physician Seating build order number form

81 ENOCHS REFERENCE

- 82 ENOCHS Warranty
- 83 ENOCHS Sales policy

84 ENOCHS CONTACT INFORMATION



ENOCHS

ENOCHS Procedural Tables

Resolute 6250 Resolute 6220 Resolute 6210



ENOCHS RESOLUTE

compare.

	RESOLUTE 6250	RESOLUTE 6220	RESOLUTE 6210				
Warranty:	4-Y	EAR	2-YEAR				
Base construction:		Fully welded all steel base					
Power System:		Electro-mechanical					
Programmable:	Yes	Programmable with optional hand control	Not available on this model				
Memory Positions:	Four	Four with optional hand control	Not available on this model				
Attainable Positions:	Chair, Contou	ured Chair, Trendelenburg, OB/GYN, Lithotomy,	Urology, Flat				
Power Adjustments:		Height: 24" to 40" Back: horizontal to 80° Tilt: horizontal to 30° Legrest: vertical to horizontal to 30°					
Auto-Return	Attainable as a programmed position	Attainable as a programmed position Attainable as a programmed position with optional hand control Not available on this mo					
Foot Control:	Easy touch, "hands-free" foot control Safety Extra-Low Voltage (SELV) and IPX1						
Hand Control:	Standard, Programmable, Safety Extra-Low Voltage (SELV)	Optional accessory that may be added in the field	Not available on this model				
Hand Control Connection Port:		oth sides of table's moveable column le, eliminating over-stretched cord	Not available on this model				
Upholstered couch Top:	One-piece style, no cross seams	One-piece style, no cross seams, embossing, or folded fabric. Securely fastened to all steel frame - no Velcro™					
Headend Accessory Rails:	Staiı	Stainless steel rails, mounted on both sides of the table					
Irrigation Basin:	2" deep stainless steel basin, inte	egrated safety switch prohibits operation of foo	t section when basin is extended				
Safety Lockout:		Accessible from both sides of table, shuts off power to all articulating sections of table, but retains power to receptacles Not available on this model					
Electrical Receptacles:	Four hospital grade electri	ical receptacles (2 per side)	Not available on this model				
Receptacle Overload Protection:	Dedicated circuit breaker	with quick reset - no fuses	Not available on this model				
Stirrups:	Positive lock, multi-pos	sition stirrups, length and lateral adjustment, co	onceal when not in use				
Paper Roll Holder:		Accommodates 18" or 21" paper widths					
Accessories:	Knee crutches, chair arms, hand control (for to footend access	6220 only), padded arm board, restraint straps, sory rails, caster base, universal clamp, multi-po	IV pole, urology basin with screen and drain, ositional clamp				
Mains Power Entry:	Master ON/OFF, IEC connection	on, Overload Protection: Integrated circuit brea	ker with quick reset - no fuses				
Duty Cycle:		1 minute ON / 9 minutes OFF					
Mode of Operation:		ttent operation tinuous operation	Table: Intermittent operation Receptacles: Not applicable				
Ratings:	Table: 115V 60 Hz 10 amps Receptacles: 115V 60 Hz Total receptacle load: 3.9 amps Table: 115V 60 Hz 10 amps Receptacles: Not applicable Total receptacle load: Not applicable						
Table Dimensions:		Length: 71", extends to 83" Width top: 30", seat tapers to 25" at headend Headrest: 19" wide, tapers to 15" Legrest: 19" wide, tapers to 18"					
Table Weight:		485 lbs.					
Patient Load:		500 lbs. evenly distributed on main couch					
Certifications:	Tested to	0 UL60601-1:2003 R6.03 CAN/CSA-C22.2 No. 60	1.1-M90				

ENOCHS



Exceptional Features:

- → 500 LBS.

 PATIENT LOAD

 evenly distributed on main couch
- → Hand Control with 4 Programmable Memory Positions
- → Power Adjustment of Height, Back, Tilt and Foot
- → Safety Lockout and 4 Electrical Receptacles

Power 6250 Table Also Features:

- 4-YEAR WARRANTY
- Safety Lockout, accessible from either side of table, immediately disables hand and foot controls to prevent accidental movement of table during exam or procedure, but maintains power to electrical receptacles
- Electrical Receptacles, 2 per side. Overload Protection: quick reset Circuit Breaker no fuses
- Receptacles, Safety Lockout, and Hand Control Ports are located on each side of table's moveable column for convenience and easy access
- Items plugged into side of table travel with the height adjustment of the table so over-stretched cords are eliminated
- "Easy Touch" Foot Control provides hands free adjustment of table's height, back, tilt and foot
- 30" Wide Seat Section provides patient more comfort
- One-piece Upholstered Top designed for easy efficient cleaning
- Adjustable Padded Headrest
- Multi-Position Positive Lock Stirrups with length and lateral adjustment
- > 2" deep Stainless Steel Irrigation Basin conveniently stores under seat section
- Stainless Steel Head-End Accessory Rails are a standard feature

Pictured above: Cream White Base, Windsor Wine Top

ENOCHS RESOLUTE 6250

Programmable Procedural Table

PATIENT LOAD: 500 LBS.

evenly distributed on the table's main couch



TABLE FEATURES:

- 4-YEAR WARRANTY
- Programmable Hand Control with 4 memory positions
- "Easy Touch" foot control for hands free table positioning
- . Power Adjustment of Back, Tilt, Foot and Height
- · Safety Lockout to prevent accidental movement
- · Circuit Breakers for overload protection-no fuses
- Electrical Receptacles, Safety Lockouts and Hand Control Ports located on table's movable column for convenience and easy access
- Cords plugged into table's receptacles, move up and down with the table, eliminating over-stretched cords
- Master ON/OFF
- Multi-Position Positive Lock Stirrups
- · Stainless Steel Irrigation Basin
- · Stainless Steel Headend Accessory Rails



Positioning includes chair, contoured chair, Trendelenburg, Lithomy, OBGYN and flat. A safety anti-pinch system stops movement of table if a pinch situation exists between headrest and floor.



Power adjustment of height (24"-40"), back (flat to 80"), tilt (flat to 30") and foot sections (vertical to horizontal) allows the 6250 to be positioned for numerous exams and procedures.

Specifications for the RESOLUTE 6250

Warranty:	4-YEAR
Base construction:	Fully welded all steel base
Power System:	Electro-mechanical
Attainable Positions:	Chair, Contoured Chair, Trendelenburg, OB/GYN, Lithotomy, Urology, Flat
Power Adjustments:	Height: 24" to 40", Back: horizontal to 80° Tilt: horizontal to 30°, Leg: vertical to horizontal
Programmable:	Quick, easy set-up of 4 memory positions
Hand Control:	Easy to operate, simple to program table positions
Hand Control Connection Port:	Located on both sides of table's moveable column for easy access and elimination of over-stretched cords
Foot Control:	"Easy touch" hands-free foot control
Upholstered Couch Top:	One-piece style, no cross seams, embossing or folded fabric, securely fastened to all steel frame
Headend Accessory Rails:	Stainless steel rails, mounted on both sides of the table
Irrigation Basin:	2" deep stainless steel basin
Electrical Receptacles:	Four hospital grade electrical receptacles (2 per side) Overload Protection: Dedicated circuit breaker with quick reset
Safety Lockout:	Shuts power off to all articulating sections of table, but retains power to receptacles
Stirrups:	Positive lock, multi-position stirrups, length and lateral adjustment
Paper Roll Holder:	Accommodates 18" or 21" paper widths
Accessories:	Knee crutches, padded arm board, restraint straps, IV pole, urology basin with screen and drain, footend accessory rails, caster base, universal clamp, multi-positional clamp
Mains Power Entry:	Master ON/OFF, IEC connection, Overload Protection: Integrated circuit breaker with quick reset
Duty Cycle:	1 minute ON / 9 minutes OFF
Ratings:	Table: 115V 60 Hz 10 amps Receptacles: 115V 60 Hz Total Receptacle Load: 3.9 amps
Table Dimensions:	Length: 71", extends to 83" with headrest Width (top): 30" at seat, tapers to 25" at headend Headrest: 19" wide, tapers to 15" Legrest: 19" wide, tapers to 18"
Table Weight:	485 lbs.
Patient Load:	500 lbs. evenly distributed on main couch
Certifications:	Tested to UL60601-1:2003 R6.03 CAN/CSA-C22.2 No. 601.1-M90







ENOCHS



Quality Features:

- → 500 LBS.
 PATIENT LOAD
 evenly distributed on main couch
- → Safety Lockout
- → Electrical Receptacles- 2 per side
- → Easy Touch
 Foot Control
- → Hand Control option available, add at any time

Resolute 6220 Table Also Features:

- ➤ 4-YEAR WARRANTY
- Electrical Receptacles, 2 per side. Overload Protection: quick reset Circuit Breaker no fuses
- Hand Control Connection Port is standard on both sides of table to allow hand control accessory to be added at any time
- Safety Lockout, accessible from either side of table, immediately disables hand and foot control to prevent accidental movement
 of table during exam or procedure, but maintains power to electrical receptacles
- Receptacles, Safety Lockout and Hand Control Ports are located on each side of table's moveable column for convenience and easy access
- Items plugged into side of table travel with the height of the table so over-stretched cords are eliminated
- "Easy Touch" Foot Control provides hands free adjustment of table's height, back, tilt and foot
- One-piece Upholstered Top designed for easy efficient cleaning.
- Adjustable Padded Headrest
- Multi-Position Positive Lock Stirrups with length and lateral adjustment
- 2" deep Stainless Steel Irrigation Basin conveniently stores under seat section.
- ➤ Stainless Steel Head-End Accessory Rails are a standard feature

Pictured above: Cream White Base, Ninja Top

ENOCHS RESOLUTE 6220

General Procedural Table

PATIENT LOAD: 500 LBS.

evenly distributed on the table's main couch



TABLE FEATURES:

4-YEAR WARRANTY

- · "Easy Touch" Foot Control for hands free adjustment of table
- Electrical Receptacles, 2 per side. Overload Protection: quick reset
 Circuit Breaker no fuses
- Safety Lockout, accessible from either side of table, prevents accidental movement of table during exam or procedure, but maintains power to electrical receptacles
- Hand Control Connection Port, standard on both sides of table, allows hand control accessory to be added at any time
- Receptacles, Safety Lockout and Hand Control Ports located on table's moveable column for convenience and easy access
- Items plugged into side of table travel with the height of the table, eliminating over-stretched cords
- · Adjustable Padded Headrest
- Multi-Position Positive Lock Stirrups
- 2" deep Stainless Steel Irrigation Basin
- Stainless Steel Headend Accessory Rails



Power adjustment of height (24"–40"), back (flat – 80"), tilt (flat – 30") and foot sections (vertical – horizontal) allows table to be positioned for various exams and procedures including chair, contoured chair, Trendelenburg, Lithorny, OBGYN and flat.



Quick access to electrical receptacles and the safety lockout switch on the side of the table saves time. The innovative Lockout System disables the hand and foot control to eliminate accidental table movement during exams or procedures but retains power to table's receptacles.

Specifications for the RESOLUTE 6220

Warranty:	4-YEAR
Base construction:	Fully welded all steel base
Power System:	Electro-mechanical
Attainable Positions:	Chair, Contoured Chair, or OB/GYN, Lithotomy, Urology, Flat
Power Adjustments:	Height: 24" to 40", Back: horizontal to 80° Tilt: horizontal to 30°, Leg: vertical to horizontal
Programmable:	Programmable with optional hand control accessory
Hand Control:	Optional accessory that may be added later
Hand Control Connection Port:	Standard feature, located on both sides of table's moveable column for easy access and elimination of over-stretched cords.
Foot Control:	"Easy touch" hands-free foot control Safety Extra-Low Voltage (SELV) and IPX1
Upholstered Couch Top:	One-piece style, no cross seams, embossing or folded fabric, securely fastened to all steel frame
Headend Accessory Rails:	Stainless steel rails, mounted on both sides of the table
Irrigation Basin:	2" deep stainless steel basin, integrated safety switch prohibits operation of foot section when basin is extended
Electrical Receptacles:	Four hospital grade electrical receptacles (2 per side) Overload Protection: Dedicated circuit breaker with quick reset
Safety Lockout:	Shuts power off to all articulating sections of table, but retains power to receptacles
Stirrups:	Positive lock, multi-position stirrups, length and lateral adjustment
Paper Roll Holder:	Accommodates 18" or 21" paper widths
Accessories:	Knee crutches, padded arm board, restraint straps, IV pole, urology basin with screen and drain, footend accessory rails, caster base, universal clamp, multi-positional clamp
Mains Power Entry:	Master ON/OFF, IEC connection, Overload Protection: Integrated circuit breaker with quick reset
Duty Cycle:	1 minute ON / 9 minutes OFF
Ratings:	Table: 115V 60 Hz 10 amps Receptacles: 115V 60 Hz Total Receptacle Load: 3.9 amps
Table Dimensions:	Length: 71", extends with headrest to 83" Width top: 30" at seat, tapers to 25" at headend Headrest: 19" wide, tapers to 15" Legrest: 19" wide, tapers to 18"
Table Weight:	485 lbs.
Patient Load:	500 lbs. evenly distributed on main couch
Certifications:	Tested to UL60601-1:2003 R6.03 CAN/CSA-C22.2 No. 601.1-M90









Distinctive Features:

- → 500 LBS.

 PATIENT LOAD

 evenly distributed on main couch
- → High Performance Height, Back, Tilt and Foot Power Adjustment
- → 30" Wide Seat Section
- → Foot Control
 - ➤ "Hands Free"
 - Utilizes Natural Left to Right Foot Movement

Resolute 6210 Table Also Features:

- 2-YEAR WARRANTY
- "Easy Touch" Foot Control provides hands free adjustment of table's height, back, tilt and foot
- One-piece Upholstered Top designed without cross seams, decorative embossing or material folds that hinder inflectional control procedures
- ➤ 30" Wide Seat Section to provide larger patients more comfort
- Mains ON/OFF shuts off all power to table to prevent accidental movement of table during exam or procedure
- ➤ Mains Power Overload Protection: easy reset Circuit Breaker no fuses.
- Adjustable Padded Headrest
- Multi-Position Positive Lock Stirrups with length and lateral adjustment.
- 2" deep Stainless Steel Irrigation Basin conveniently stores under seat section.
- ➤ Stainless Steel Head-End Accessory Rails are a standard feature

ENOCHS RESOLUTE 6210

Basic Procedural Table

PATIENT LOAD: 500 LBS.

evenly distributed on the table's main couch



TABLE FEATURES:

• 2-YEAR WARRANTY

- "Easy Touch" Foot Control for hands free adjustment of table
- One-piece Upholstered Top designed without cross seams, decorative embossing or material folds that hinder inflectional control procedures
- 30" Wide Seat Section to provide larger patients more comfort
- Mains ON/OFF shuts off all power to table to prevent accidental movement of table during exam or procedure
- Mains Power Overload Protection: easy reset Circuit Breaker - no fuses
- Adjustable Padded Headrest
- Multi-Position Positive Lock Stirrups
- 2" deep Stainless Steel Irrigation Basin
- Stainless Steel Headend Accessory Rails



Positive lock, multi-position stirrups with lateral and length adjustment and a 2" deep stainless steel irrigation basin that quickly and easily removes for cleaning after each procedure are standard features of the 6210. Table is engineered with an integrated safety switch that stops the leg section from going UP when irrigation basin is not in full stored position, eliminating damage to the basin and/or leg extension.



The table's 30" wide seat section provides maximum comfort to larger patients while the one-piece seamless top surface, designed without decorative embossing or material creases and folds that collect dirt, allows staff to quickly and efficiently clean after each patient. The fully padded upholstered top, covered in durable Naugahyde®, is securely attached to an all steel frame for strength and stability.

Specifications for the RESOLUTE 6210

Warranty:	2-YEAR
Base construction:	Fully welded all steel base
Power System:	Electro-mechanical
Attainable Positions:	Chair, Contoured Chair, Trendelenburg, OB/GYN, Lithotomy, Urology, Flat
Power Adjustments:	Height: 24" to 40", Back: horizontal to 80° Tilt: horizontal to 30°, Leg: vertical to horizontal
Programmable:	Not available on this model
Hand Control:	Not available on this model
Hand Control Connection Port:	Not available on this model
Foot Control:	"Easy touch" hands-free foot control Safety Extra-Low Voltage (SELV) and IPX1
Upholstered Couch Top:	One-piece style, no cross seams, embossing or folded fabric, securely fastened to all steel frame
Headend Accessory Rails:	Stainless steel rails, mounted on both sides of the table
Irrigation Basin:	2" deep stainless steel basin, integrated safety switch prohibits operation of foot section when basin is extended
Electrical Receptacles:	Not available on this model
Safety Lockout:	Not available on this model
Stirrups:	Positive lock, multi-position stirrups, length and lateral adjustment
Paper Roll Holder:	Accommodates 18" or 21" paper widths
Accessories:	Knee crutches, padded arm board, restraint straps, IV pole, urology basin with screen and drain, footend accessory rails, caster base, universal clamp, multi-positional clamp
Mains Power Entry:	Master ON/OFF, IEC connection, Overload Protection: Integrated circuit breaker with quick reset
Duty Cycle:	1 minute ON / 9 minutes OFF
Ratings:	Table: 115V 60 Hz 10 amps Receptacles: N/A , Total receptacle load: N/A
Table Dimensions:	Length: 71", extends with headrest to 83" Width top: 30" at seat, tapers to 25" at headend Headrest: 19" wide, tapers to 15" Legrest: 19" wide, tapers to 18"
Table Weight:	485 lbs.
Patient Load:	500 lbs. evenly distributed on main couch
Certifications:	Tested to UL60601-1:2003 R6.03 CAN/CSA-C22.2 No. 601.1-M90







RESOLUTE Power Series Accessories

Model #6109 Universal Clamp




Model #6110 Tri-Clamp



Model #6112 Base Rails (pair)



Model #6122 Knee Crutches (pair)

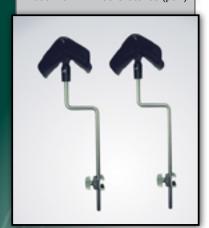




Requires #6109 or 6110 Clamp

Model #6124 Padded Armboard

Model #6126 Restraint Straps



Requires #6112 Base Rails

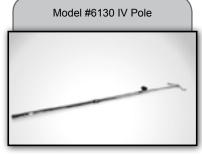

Requires #6109 or 6110 Clamp

Model #6140 Chair Arms

Model #6145 Programmable Hand Control



Model #6137 Urology with Screen, Drain & Hose



Model #6150 Caster Base





RESOLUTE Color Selections

UPHOLSTERY COLORS



Ninja



Smoke



Cambridge Blue



Simply Taupe



Windsor Wine



Newberry Navy



Ocean Jade



Fawn



Sussex Green



British Tan



India Ink

BASE COLORS




Pearl Gray



Cream White

RESOLUTE Series of Procedural Tables

Build Order Number

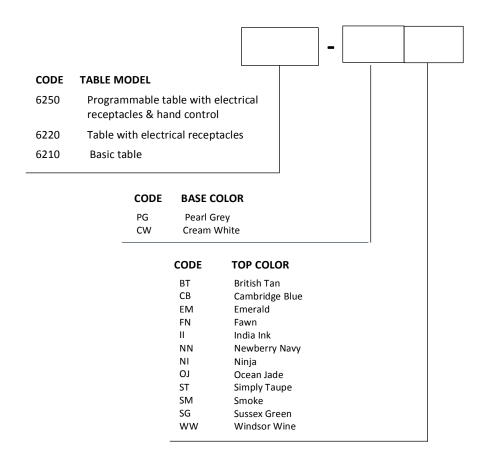

ENOCHS

ENOCHS Low-Access High-Low Tables

> Encore 4550 Encore 4520 Encore 4510



ENOCHS ENCORE

LOW ACCESS

compare.

	ENCORE 4550 Programmable Low Access Table	ENCORE 4520 Low Access Table	ENCORE 4510 Basic Low Access Table				
Warranty:	4-YEAR	4-YEAR	2-YEAR				
Operation:	Programmable with 3 memory positions	Non-programmable	Non-programmable				
Patient Load:	700 Lbs.						
Power Adjustments:	He	eight: Low-Access: 18.5" to 38" maximum Back: Horizontal to 80°					
Controls:	Hand Control, Programmable	Foot Control	Foot Control				
Couch Top:	100% seamle	ss vacuum formed style, all steel frame construction	on				
Perineal Recess:	Moul	lded, no seams, 10" W x 2" D at center point					
Pelvic Tilt:		11°					
Leg Extension:	E	xtends out 20" from end of seat section					
Irrigation Basin:		Stainless Steel, 2.5" deep					
Master ON/OFF Switch:	Shuts off power to entire product						
Stirrups:	Multi-positional, lateral, and length adjustment						
Storage:	Side Drawer: Field reve	rsible, 15.5" W x 15" D x 4" H	No side drawer				
Storage.	Footend Drawer: stainless steel, 16" W x 10.75" D x 4" H						
Base Construction:		Fully welded all steel base					
Paper Roll Holder:	Acc	ommodates 18" or 21" table paper widths					
Accessories:		nes, padded armboard, drawer locks, restraint stra n screen and drain, headend accessory rails, foot e					
Power System:		Electro-mechanical operating system					
Mains Power Entry:	IEC connec	tion with internal circuit breaker, detachable cord					
Table Duty Cycle:		1 minute run/9 minutes rest					
Table Rating:		115V 60Hz 10 amps					
Table Dimensions:	Width: 27'	", Height: 18.5" - 38", Length: 52" (extends to 72")					
Certifications:	Classified by UL to ANSI/AAMI	ES60601-1 (2005) + AMD 1 (2012), CAN/CSA-C22.2 N	lo. 60601-1 (2014)				
	ENCORE 4550 Add-On Options	ENCORE 4520 Add-On Options	ENCORE 4510 Add-On Option				

ENCORE 4550 Add-On Options					ENCORE	4520 Add-On (E	NCORE 4510 Add-On Options	
	Electrical Receptacles*	Instrument Warmer**	Foot Control		Electrical Receptacles*	Instrument Warmer**	Hand Control: non-programmable		No Add-On Options
-01	•	•		-01	•	•			
-02	•			-02	•				
-03		•		-03		•			
-04	•	•	•	-04	•	•	•		
-05		•	•	-05		•	•		
-06	•		•	-06	•		•		
-07			•	-07			•		

^{*} Hospital grade, dual receptacles with easy access dedicated circuit breaker **Duty Cycle:** Continuous operation **Total Load:** 3.9 amps Rating: 115V 60Hz

^{** 45} Watts, ON/OFF switch with integrated circuit breaker **Duty Cycle:** Continuous operation

Programmable Low-Access Examination Table





New Features:

- → 18.5" LOW ACCESS
 - ➤ ADA compliant
 - ➤ 38" height max.
- → 700 LBS.

 PATIENT LOAD
- → Programmable with 3 memory positions
- → Power Back
- → 11° Pelvic Tilt
- → Reversible Side Drawer

ENCORE Power 4550 Programmable Low-Access Table Also Features:

➤ 4-YEAR WARRANTY

- > Programmable Hand Control
- ➤ All Steel Base Construction NO PLASTIC
- ➤ 100% Seamless Vacuum Formed Top with moulded Perineal Recess
- ➤ Multi-position Stirrups with length and lateral adjustment
- ➤ Pull out padded Leg Extension
- ➤ 2" deep stainless steel irrigation basin moves independent of Leg Extension
- > Stainless Steel foot end drawer
- ➤ ADD ON OPTIONS:
 - Electrical Receptacles
 - Instrument Warming Drawer
 - Easy Touch Foot Control

ENOCHS ENCORE 4550

PROGRAMMABLE LOW-ACCESS EXAMINATION TABLE

PATIENT LOAD: 700 LBS.

evenly distributed on the table's main couch



TABLE FEATURES:

- 4-YEAR Warranty
- Patient Load: 700 lbs.
- Programmable Hand Control with 3 memory positions
- Low access: 18.5" Maximum height: 38"
- Power Back adjusts from horizontal to 80°
- Easy to clean 100% Seamless Vacuum Formed Top with additional padding on edges for patient comfort
- 11° Pelvic Tilt
- Field Reversible Side Drawer
- Stainless Steel Footend Drawer
- Stainless Steel Irrigation Basin operates independent of Leg Extension
 Multi-position Positive Lock Stirrups with length and lateral adjustment
 Overload Protection: Circuit Breakers NO FUSES
- Fully Welded All Steel Base Construction NO PLASTIC
- Classified by UL to ANSI/AAMI ES60601-1 (2005) + AMD 1 (2012), CAN/CSA-C22.2 No. 60601-1 (2014)



Working end of the table features multiposition stirrups and a padded pullout leg extension designed to move independent of the stainless steel irrigation basin for increased efficiency. The stainless steel footend drawer provides storage and accommodates the instrument warmer when ordered as an ADD-ON OPTION.



The ENOCHS 4550 is ADA compliant with safe, easy, no barrier access onto table's 18.5" low seat height. The power back adjusts to any angle from horizontal to 80° allowing couch surface to be flat or in a chair position for confident patient transfer. A Grab Rail System for patient assistance is available as an accessory.

Specifications for the ENCORE 4550 Programmable Low-Access Table

Warranty:	4-YEAR			
Operation:	Power, programmable with 3 memory positions			
Patient Load:	700 lbs. evenly distributed on main couch			
Power Adjustments:	Height: Low Access: $18.5'' \pm .25''$ to maximum $38'' \pm .25''$ Back: Horizontal to $80^{\circ} \pm 2^{\circ}$			
Controls:	"Easy Read" hand control, programmable, SELV			
Couch Top:	100% seamless vacuum formed style, all steel frame construction			
Perineal Recess:	Moulded, no seams, 10" W x 2" D at center point			
Pelvic Tilt:	11° ± 2°			
Leg Extension:	Extends out 20" from end of seat section			
Irrigation Basin:	Stainless steel , 2.5" deep			
Master ON/OFF Switch:	Shuts off power to entire product			
Stirrups:	Multi-positional, lateral, and length adjustment			
Storage:	Side Drawer: Field reversible, 15.5" W x 15" D x 4" H Footend Drawer: Stainless Steel, 16" W x 10.75" D x 4" H			
Base Construction:	Fully welded all steel base			
Paper Roll Holder:	Accommodates 18" or 21" table paper widths			
Accessories:	Grab rails, knee crutches, padded armboard, drawer lock, restraint straps, IV pole, urology basin with screen and drain, headend accessory rails, foot end rails			
Power System:	Electro-mechanical operating system			
Mains Power Entry:	IEC connection with internal circuit breaker, detachable cord			
Table Duty Cycle:	1 minute run / 9 minutes rest			
Table Rating:	115V 60Hz 10 amps			
Table Dimensions:	Width: 27", Height: 18.5" - 38", Length: 52" (extends to 72")			
Certifications:	Classified by UL to ANSI/AAMI ES60601-1 (2005) + AMD 1 (2012), CAN/CSA-C22.2 No. 60601-1 (2014)			

	AD	D-ON OPTIONS	
	Electrical Receptacles	Instrument Warmer	Foot Control
-01	•	•	
-02	•		
-03		•	
-04	•	•	•
-05		•	•
-06	•		•
-07			•
	Hospital grade, dual receptacles with easy access dedicated circuit breaker Duty Cycle: Continuous operation Total Load: 3.9 amps Rating: 115V 60Hz	45 Watts, 0N/OFF switch with integrated circuit breaker, heats to approximate body temperature Duty Cycle: Continuous operation	"Hands Free", easy to operate, SELV, non-detachab







ENCORE 4550 Programmable Low Access Exam Table

Build Order Number

				4	550	- [
CODE		ADD-ON OPTIO	NC.									
00		ble specified in										
01				nd Instrument W	armer							
02		al Electrical Red										
03	add Inst	trument Warm	er only	•								
04	add Dua	al Electrical Red	ceptacles, I	nstrument Warm	ner and Foot	Con	trol					
05		trument Warm										
06			ceptacles a	nd Foot Control								
07	add Foo	ot Control only										
	CODE	SIDE DRAWER		TION				•				
	RH	Patient's RIG										
	LH	Patient's LEF	T side									
		CODE	DRAWER	FRONT COLOR								
		GR	Grey									
		BE	Beige									
			CODE	BASE COLOR								
			PG	Pearl Grey								
			CW	Cream White								
				CODE	TOP COLO	OR						
				СВ	Cambridge	e Blue	9					
				EB	English Bro							
				HS	Highland S	age						
				II	India Ink							
				NN	Newberry		′					
				OJ	Ocean Jad							
				ST	Simply Tau	ıpe						
				SM	Smoke							
				WW	Windsor V	Vine						

Low-Access Examination Table





New Features:

- → 18.5" LOW ACCESS
 - ➤ ADA compliant
 - > 38" height max.
- → 700 LBS.

 PATIENT LOAD
- → Power Back
- → "Easy Touch" Foot Control
- → 11° Pelvic Tilt
- → Reversible Side Drawer

ENCORE Power 4520 Low-Access Table Also Features:

➤ 4-YEAR WARRANTY

- ➤ Easy Touch Foot Control
- ➤ All Steel Base Construction NO PLASTIC
- ➤ 100% Seamless Vacuum Formed Top with moulded Perineal Recess
- ➤ Multi-position Stirrups with length and lateral adjustment
- ➤ Pull out padded Leg Extension
- ➤ 2" deep stainless steel irrigation basin moves independent of Leg Extension
- > Stainless Steel foot end drawer
- ➤ ADD ON OPTIONS:
 - Electrical Receptacles
 - Instrument Warming Drawer
 - Standard Hand Control (non-programmable)

ENOCHS ENCORE 4520

LOW-ACCESS EXAMINATION TABLE

PATIENT LOAD: 700 LBS.



TABLE FEATURES:

- 4-YEAR Warranty
- Patient Load: 700 lbs.
- Low access: 18.5" Maximum height: 38"
- Power Back adjusts from horizontal to 80°
- Easy to clean 100% Seamless Vacuum Formed Top with additional padding on edges for patient comfort
- 11° Pelvic Tilt
- Field Reversible Side Drawer
- Stainless Steel Footend Drawer
- Stainless Steel Irrigation Basin operates independent of Leg Extension
- Multi-position Positive Lock Stirrups with length and lateral adjustment
- Overload Protection: Circuit Breakers NO FUSES
- Fully Welded All Steel Base Construction NO PLASTIC
- Classified by UL to ANSI/AAMI ES60601-1 (2005) + AMD 1 (2012), CAN/CSA-C22.2 No. 60601-1 (2014)



Facilities create their perfect table by selecting ADD-ON OPTIONS they want, not forced to buy features they do not need. Options of hospital grade dual electrical receptacles with dedicated circuit breaker, instrument warmer with easy access ON/OFF switch and a standard non-programmable SELV hand control are all available to order with a table.



The ENOCHS 4520 is ADA compliant offering patients easy, no barrier access to table. The power back adjusts to any angle from horizontal to 80° and allows couch surface to be flat or in a chair position for patient transfers. A Grab Rail System to provide assistance for the patient is available as an accessory. The "Easy Touch" SELV foot control comes with the 4520 allowing hands free adjustment of table. A non-programmable hand control may be chosen as an ADD-ON OPTION or purchased at a later date.

Specifications for the ENCORE 4520 General Low-Access Table

Warranty:	4-YEAR			
Operation:	Power, non-programmable			
-				
Patient Load:	700 lbs. evenly distributed on main couch			
Power Adjustments:	Height: Low Access: $18.5'' \pm .25''$ to maximum $38'' \pm .25''$ Back: Horizontal to $80^{\circ} \pm 2^{\circ}$			
Controls:	"Easy Touch" hands free, Foot Control, SELV			
Couch Top:	100% seamless vacuum formed style, all steel frame construction			
Perineal Recess:	Moulded, no seams, 10" W x 2" D at center point			
Pelvic Tilt:	11° ± 2°			
Leg Extension:	Extends out 20" from end of seat section			
Irrigation Basin:	Stainless steel , 2.5" deep			
Master ON/OFF Switch:	Shuts off power to entire product			
Stirrups:	Multi-positional, lateral, and length adjustment			
Storage:	Side Drawer: Field reversible, 15.5" W x 15" D x 4" H Footend Drawer: Stainless Steel, 16" W x 10.75" D x 4" H			
Base Construction:	Fully welded all steel base			
Paper Roll Holder:	Accommodates 18" or 21" table paper widths			
Accessories:	Grab rails, knee crutches, padded armboard, drawer lock, restraint straps, IV pole, urology basin with screen and drain, headend accessory rails, foot end rails			
Power System:	Electro-mechanical operating system			
Mains Power Entry:	IEC connection with internal circuit breaker, detachable cord			
Table Duty Cycle:	1 minute run / 9 minutes rest			
Table Rating:	115V 60Hz 10 amps			
Table Dimensions:	Width: 27", Height: 18.5" - 38", Length: 52" (extends to 72")			
Certifications:	Classified by UL to ANSI/AAMI ES60601-1 (2005) + AMD 1 (2012), CAN/CSA-C22.2 No. 60601-1 (2014)			

	AD	D-ON OPTIONS	
	Electrical Receptacles	Instrument Warmer	Hand Control
-01	•	•	
-02	•		
-03		•	
-04	•	•	•
-05		•	•
-06	•		•
-07			•
	Hospital grade, dual receptacles with easy access dedicated circuit breaker Duty Cycle: Continuous operation Total Load: 3.9 amps Rating: 115V 60Hz	45 Watts, ON/OFF switch with integrated circuit breaker, heats to approximate body temperature Duty Cycle: Continuous operation	"Easy Read", non-programmable, SELV








enochs@enochsmed.com

ENCORE 4520 Low Access Exam Table

Build Order Number

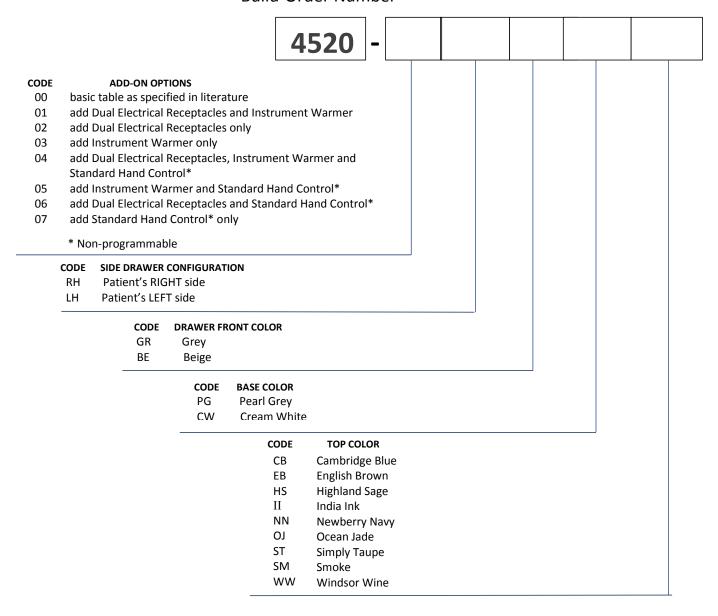

Basic Low-Access Examination Table





New Features:

- → 18.5" LOW ACCESS
 - ➤ ADA compliant
 - > 38" height max.
- → 700 LBS.

 PATIENT LOAD
- → Power Back
- → "Easy Touch" Foot Control
- → 11° Pelvic Tilt

ENCORE Power 4510 Basic Low-Access Table Also Features:

➤ 2-YEAR WARRANTY

- ➤ Easy Touch Foot Control
- ➤ All Steel Base Construction NO PLASTIC
- ➤ 100% Seamless Vacuum Formed Top with moulded Perineal Recess
- ➤ Multi-position Stirrups with length and lateral adjustment
- > Pull out padded Leg Extension
- ➤ 2" deep irrigation basin moves independent of Leg Extension
- > Stainless Steel foot end drawer
- ➤ No ADD ON OPTIONS are available on this model

ENOCHS ENCORE 4510

BASIC LOW-ACCESS EXAMINATION TABLE

PATIENT LOAD: 700 LBS.



TABLE FEATURES:

- 2-YEAR Warranty
- Patient Load: 700 lbs.
- Low access: 18.5" Maximum height: 38"
- Power Back adjusts from horizontal to 80°
- Easy to clean 100% Seamless Vacuum Formed Top with additional padding on edges for patient comfort
- 11° Pelvic Til
- Stainless Steel Footend Drawer
- Stainless Steel Irrigation Basin operates independent of Leg Extension
- Multi-position Positive Lock Stirrups with length and lateral adjustment
- Overload Protection: Circuit Breakers NO FUSES
- Fully Welded All Steel Base Construction NO PLASTIC
- Classified by UL to ANSI/AAMI ES60601-1 (2005) + AMD 1 (2012), CAN/CSA-C22.2 No. 60601-1 (2014)



back with reliable hands free action. Natural left to right rocking motion of operator's foot makes table adjustments effortless.

Multi-position stirrups, a padded pull-out leg extension and a separate pull-out Stainless Steel irrigation basin are all features of the 4510 designed to reduce time and increase efficiency in exams or procedures.



The 2" deep recess, moulded into the seat section, provides more comfortable positioning during perineal exams or procedures. The 11° pelvic tilt with release handles located on both sides of the table is a standard feature.



A Grab Rail System accessory for ENCORE tables provides assistance to patients when getting on or off the table. The unique design of ENOCHS' all steel grab rails allows staff to quickly remove one or both units when not needed and easily reinstall when required.

Specifications for the ENCORE 4510 Basic Low-Access Table

Warranty:	2-YEAR
Operation:	Power, non-programmable
Patient Load:	700 lbs. evenly distributed on main couch
Power Adjustments:	Height: Low Access: $18.5'' \pm .25''$ to maximum $38'' \pm .25''$ Back: Horizontal to $80^{\circ} \pm 2^{\circ}$
Controls:	"Easy Touch" hands free, Foot Control, SELV
Couch Top:	100% seamless vacuum formed style, all steel frame construction
Perineal Recess:	Moulded, no seams, 10" W x 2" D at center point
Pelvic Tilt:	11° ± 2°
Leg Extension:	Extends out 20" from end of seat section
Irrigation Basin:	Stainless steel , 2.5" deep
Master ON/OFF Switch:	Shuts off power to entire product
Stirrups:	Multi-positional, lateral, and length adjustment
Storage:	Side Drawer: Not available on this model Footend Drawer: Stainless Steel, $16''$ W x $10.75''$ D x $4''$ H
Base Construction:	Fully welded all steel base construction
Paper Roll Holder:	Accommodates 18" or 21" table paper widths
Power System:	Electro-mechanical operating system
Accessories:	Grab rails, knee crutches, padded armboard, restraint straps, IV pole, urology basin with screen and drain, headend accessory rails, foot end rails
Mains Power Entry:	IEC connection with internal circuit breaker, detachable cord
Table Duty Cycle:	1 minute run / 9 minutes rest
Table Rating:	115V 60Hz 10 amps
Table Dimensions:	Width: 27", Height: 18.5" - 38", Length: 52" (extends to 72")
Certifications:	Classified by UL to ANSI/AAMI ES60601-1 (2005) + AMD 1 (2012), CAN/CSA-C22.2 No. 60601-1 (2014)

No Add-on Options Available





ENCORE 4510 Basic Low Access Exam Table Build Order Number

			4510 -		ı	
CODE	DRAWER I	FRONT COLOR				
GR	Grey					
BE	Beige					
	CODE	BASE COLOR				
	PG	Pearl Grey				
	CW	Cream White				
		CODE	TOP COLOR		_	
		СВ	Cambridge Blue			
		EB	English Brown			
		HS	Highland Sage			
		II	India Ink			
		NN	Newberry Navy			
		OJ	Ocean Jade			
		ST	Simply Taupe			
		SM	Smoke			
		WW	Windsor Wine			



ENCORE 4500 Series Accessories



With Adaptor

Non-programmable. Available on 4520 only.

ENOCHS

ENOCHS Easy-Access High-Low Tables

> Encore 4250 Encore 4200



ENOCHS ENCORE

EASY ACCESS

compare.

	ENCORE 4250	ENCORE 4200				
Warranty:	4 Year	2 Years				
Base Construction:	Fully welded all steel base					
Power System:	Electro-mechanical operating system					
Power Adjustments:	Height: 24" – 40"					
•	Back: Horizontal to 80°					
Foot Control:	Hands free, standard on every table, SELV					
Hand Control:	Easy to read, optional accessory, add at any time, SELV					
Couch Top:	100% seamless vacuum formed style, all steel frame construction					
Perineal Recess:	Moulded, no seams, 10" W x 2" D at center point					
Pelvic Tilt:	7°,	11°				
Leg Extension:	Extends out 20" from end of seat section					
Irrigation Basin:	4" Deep stainless steel basin					
Electrical Receptacles:	Hospital grade, UL Listed, dual electrical outlets with dedicated circuit breaker	Not available on this model				
Auto Return/Stop:	Not available on this model					
Instrument Warmer: with ON/OFF Switch, Circuit Breaker, & Reset	High efficiency warmer utilizing footend drawer	Not available on this model				
Master ON/OFF Switch:	Shuts off power to entire product					
Stirrups:	irrups: Multi-positional, lateral, and length adjustment					
Storage:	Side Drawer: 15.5" W x 15" D x 6" H Footend Drawer: 15.5" W x 15" D x 4" H Total Storage Space: 1.4 cubic feet					
Paper Roll Holder:	Accommodates 18" or 21" table paper widths					
Accessories:	hand control, knee crutches, padded arm board, drawer locks, restraint straps, IV pole, urology basin with screen					
Mains Power Entry:	IEC connection with internal circuit breaker, detachable cord					
Duty Cycle:	Table: 1 minute ON/9 minutes OFF Receptacles: Continuous operation Instrument Warming Drawer: Continuous operation	Table: 1 minute ON/9 minutes OFF Receptacles: Not applicable Instrument Warming Drawer: Not applicable				
Rating:	Table: 115V 60Hz 10 amps Receptacles: 115V 60Hz Total receptacle load: 3.9 amps	Table: 115V 60Hz 10 amps Receptacles: Not applicable Total receptacle load: Not applicable				
Table Dimensions:	Table Dimensions: 27" Wide x 18.5"-38" High x 52.5" Long (extends to 72.5")					
Table Weight:	Table Weight: 416 Lbs.					
Patient Load:	500 Lbs., evenly distributed on main couch					
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1					

ENOCHS ENCORE POWER 4250

High-Low Examination Table

ENOCHS EXAMINING ROOM FURNITURE



Encore Power 4250 Table Also Features:

- 4-YEAR WARRANTY
- Use of Circuit Breakers for easy reset NO FUSES
- All Steel Base Construction NO PLASTIC
- Hand Control available as optional accessory
- Foot Control, standard on every table, Safety Extra Low Voltage (SELV)
- 100% Seamless Vacuum Formed Top with moulded Perineal Recess
- Pelvic Tilt 7° and 11°
- Multi-Position Stirrups with length and lateral adjustment
- Pull out padded Leg Extension
- 4" deep Stainless Steel Irrigation Basin moves independent of Leg Extension
- 2 Storage Drawers with a total of 1.4 cubic feet of storage

Pictured above: Grey Fronts, Pearl Grey Base, Smoke Top

New Features:

- → 500 LBS.

 PATIENT LOAD
- → High Performance Actuator System
 - Power Height
 - Power Back
 - ➤ Smooth & Quiet
- → Field Reversible Side Drawer
- → Electrical Outlets
 - ➤ Dual, Hospital Grade
 - Dedicated Circuit Breaker
- → Instrument Warmer
 - Efficient Warming Technology
 - ON/OFF Switch with Integrated Circuit Breaker

ENOCHS ENCORE 4250

High-Low Examination Table

PATIENT LOAD: 500 LBS.

evenly distributed on the table's main couch



'EASY TOUCH' FOOT CONTROL adjusts table's height and back with reliable hands-free action. Pedals are covered with non-slip tread and utilize a natural left to right rocking motion of operator's foot.

TABLE FEATURES:

- 4-YEAR Warranty
- Hand Control available as an accessory
- Fully welded All Steel Base Construction -NO PLASTIC
- Full Power adjustment of back
- Easy to clean 100% Seamless Vacuum Formed Top with additional padding on edges for patient comfort
- Pelvic Tilt 7° and 11°
- Electrical Outlets with circuit breaker NO FUSES
- Instrument Warming Drawer efficiently warms drawer contents to approximate body temperature
- Instrument Warmer ON/OFF Switch with Integrated Circuit Breaker - NO FUSES
- Stainless Steel Irrigation Basin operates independent of Leg Extension
- Easy to adjust Stirrups with length and lateral adjustment
- 1.4 cubic feet of Storage



The hospital grade dual electrical outlets are conveniently located at the foot end of the table. The outlets feature a dedicated circuit breaker that will trip if outlets are overloaded. A quick reset of the circuit breaker re-establishes power to the outlets, eliminating costly down time. The Power 4250 also features a SELV (safety extra low voltage) hand control, an accessory that can be ordered at time of purchase or added at a later date.



The footend of the table features a padded pullout leg extension, 4" deep stainless steel irrigation basin and the next generation instrument warming drawer system. The exclusive design of separate operation between the leg extension and basin allows for fast transition from one to other. The pull-out leg extension enables patients with limited hip rotation to be placed in a GYN position, while still giving the physician unrestricted access to the irrigation basin. The instrument warming drawer is designed to warm speculums and instruments to approximate body temperature.

Warranty:	4-YEAR		
Base construction:	Fully welded all steel base		
Power System:	Electro-Mechanical provides smooth and quiet movement		
Power Adjustments:	Height: 24" to 40" Back: horizontal to 80°		
Foot Control:	"Hands free", standard on every table, SELV		
Hand Control:	"Easy to read", optional accessory, add at any time, SELV		
Couch Top:	100% Seamless vacuum formed style, all steel frame construction		
Perineal Recess:	Moulded, no seams, 10" W x 2" D at center point		
Pelvic Tilt:	7°, 11°		
	Extends out 20" from end of seat section		
Leg Extension:			
Irrigation Basin:	4" Deep stainless steel basin		
Electrical Outlets: with Circuit Breaker	Standard feature		
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	Standard feature		
Stirrups:	Multi-positional, lateral and length adjustment		
	Side Drawer: 15.5: W x 15" D x 6" H		
Storage:	Footend Drawer: 15.5: W x 15" D x 4" H		
	Total Storage Space: 1.4 cibic feet		
Paper Roll Holder:	Accommodates 18" or 21" paper widths		
Accessories:	hand control, knee crutches, padded arm board, drawer locks, restraint straps, IV pole, urology basin with screen and drain, headend accessory rails, footend accessory rails, caster base		
Mains Power Entry: with Circuit Breaker	IEC connection, internal circuit breaker with reset		
	Table: 1 minute ON / 9 minutes OFF		
Duty Cycle:	Receptacles: Continuous operation		
	Instrument Warming Drawer: Continuous operation		
	Table: 115V 60 Hz 10 amps		
Ratings:	Receptacles: 115V 60Hz		
Takla Dimanaia ya	Total Receptacle Load: 3.9 amps		
Table Dimensions:	Width: 27", Height: 24" - 40", Length: 52" (extends to 72")		
Table Weight:	416 lbs.		
Patient Load:	500 lbs. evenly distributed on main couch		
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1		





ENOCHS ENCORE POWER 4200

High-Low Examination Table

ENOCHS EXAMINING ROOM FURNITURE



Encore Power 4200 Table Also Features:

2-YEAR WARRANTY

- Use of Circuit Breakers for easy reset NO FUSES
- All Steel Base Construction NO PLASTIC
- Hand Control available as optional accessory
- Foot Control, standard on every table, Safety Extra Low Voltage (SELV)
- 100% Seamless Vacuum Formed Top with moulded Perineal Recess
- Pelvic Tilt 7° and 11°
- Multi-Position Stirrups with length and lateral adjustment
- Pull out padded Leg Extension
- 4" deep Stainless Steel Irrigation Basin moves independent of Leg Extension
- 2 Storage Drawers with a total of 1.4 cubic feet of storage

Pictured above: Grey Fronts, Pearl Grey Base, Cambridge Blue Top

New Features:

- → 500 LBS.

 PATIENT LOAD
- → High Performance Actuator System
 - ➤ Power Height
 - > Power Back
 - ➤ Smooth & Quiet
- → Field Reversible Side Drawer

→ Foot Control


- Standard Feature
- ➤ Safety Extra Low Voltage (SELV)
- Utilizes Natural Left to Right Foot Movement

ENOCHS ENCORE 4200

High-Low Examination Table

PATIENT LOAD: 500 LBS.

evenly distributed on the table's main couch

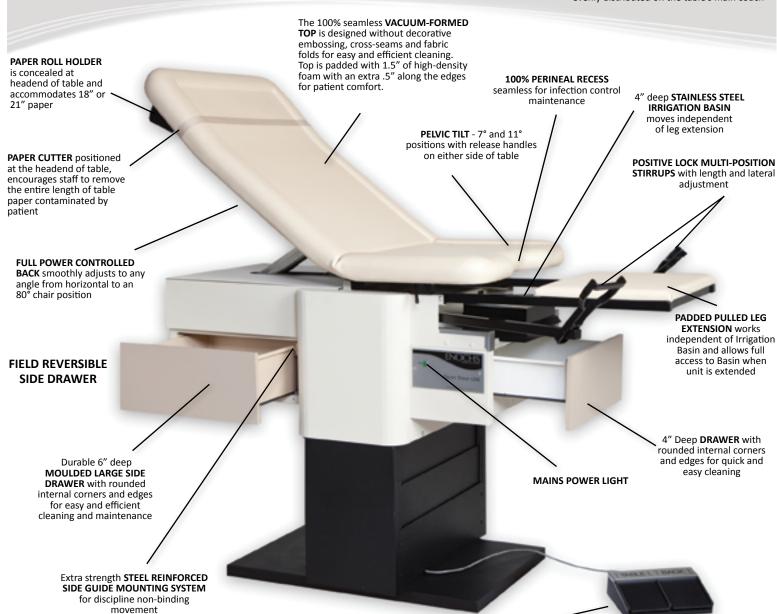

'EASY TOUCH' FOOT CONTROL adjusts table's height and back with reliable hands-free action. Pedals are covered with non-slip tread and utilize a natural left to right rocking motion of operator's foot.

TABLE FEATURES:

- 2-YEAR Warranty
- Hand Control available as an accessory
- Fully welded All Steel Base Construction NO PLASTIC
- Full Power adjustment of back
- Easy to clean 100% Seamless Vacuum Formed Top with additional padding on edges for patient comfort
- Pelvic Tilt 7° and 11°
- Stainless Steel Irrigation Basin operates independent of Leg Extension
- Easy to adjust Stirrups with length and lateral adjustment
- 1.4 cubic feet of Storage



Securely attached to the table base, the ENOCHS vacuum-formed top is the strongest top on the market with an all steel frame and back. Top is padded for patient comfort with 1.5" of high density foam with an additional .5" of foam contouring the edges. The sleek Euro-look top features a seamless surface, moulded perineal recess and a pelvic tilt with positions of 7° and 11°. The top offers proven performance of strength, durability and long service life.



The working end of the table has the exclusive ENOCHS design of the pullout padded leg extension moving independent of the 4" deep stainless steel irrigation basin. This design reduces surface area to allow fast easy clean up after irrigation procedure. Patients with limited rotation can be placed in a GYN position utilizing the pad board instead of stirrups while physician maintains unrestricted access to the irrigation basin.


Warranty:	2-YEAR	
Base construction:	Fully welded all steel base	
Power System:	Electro-Mechanical provides smooth and quiet movement	
Power Adjustments:	Height: 24" to 40"	
Foot Control:	Back: horizontal to 80° "Hands froe" standard on every table. SELV	
	"Hands free", standard on every table, SELV	
Hand Control:	"Easy to read", optional accessory, add at any time, SELV	
Couch Top:	100% Seamless vacuum formed style, all steel frame construction	
Perineal Recess:	Moulded, no seams, 10" W x 2" D at center point	
Pelvic Tilt:	7°, 11°	
Leg Extension:	Extends out 20" from end of seat section	
Irrigation Basin:	4" Deep stainless steel basin	
Electrical Outlets: with Circuit Breaker	Not available on this model	
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	Not available on this model	
Stirrups:	Multi-positional, lateral and length adjustment	
	Side Drawer: 15.5: W x 15" D x 6" H	
Storage:	Footend Drawer: 15.5: W x 15" D x 4" H	
	Total Storage Space: 1.4 cibic feet	
Paper Roll Holder:		
Accessories:	hand control, knee crutches, padded arm board, drawer locks, restraint straps, IV pole, urology basin with screen and drain, headend accessory rails, footend accessory rails, caster base	
Mains Power Entry: with Circuit Breaker	IFC connection internal circuit breaker with reset	
	Table: 1 minute ON / 9 minutes OFF	
Duty Cycle:	Receptacles: N/A	
	Instrument Warming Drawer: N/A	
	Table: 115V 60 Hz 10 amps	
Ratings:	Receptacles: N/A	
Takla Dimamaia ya	Total Receptacle Load: N/A	
Table Dimensions:	Width: 27", Height: 24" - 40", Length: 52" (extends to 72")	
Table Weight:	416 lbs.	
Patient Load:	500 lbs. evenly distributed on main couch	
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1	





ENCORE 4200 Series Easy Access Exam Table

Build Order Number Example

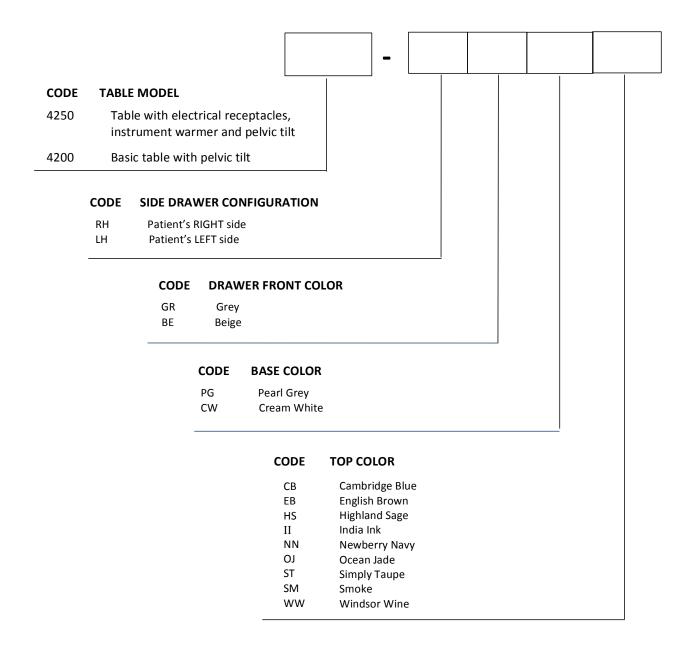



ENCORE 4200 Series Accessories


Model #6109 Universal Clamp



Model #6110 Tri-Clamp



Model #4111 Headend Rails

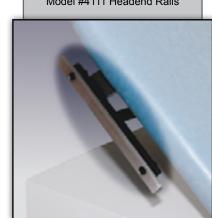

Model #4112 Footend Rails



Model #4124 Arm Board



Requires #4111 Headend Rails and #6109 or #6110 clamp









Model #4150 Caster Base



Model #4137 Urology Pan



Model #4145 Hand Control

Requires #4112 Base Rails







4520 in Smoke



4200 in Ocean Jade



Cambridge Blue



Windsor Wine



English Brown



Desert Sand



Smoke (above)



Ocean Jade (above)



Highland Sage



India Ink



Newberry Navy



Simply Taupe

BASE/DRAWER FRONT COLORS



Cream White Base Beige Drawer Fronts



Pearl Gray Base Gray Drawer Fronts



Cream White Base Natural Oak Drawer Fronts *Only available on the ENCORE 4200 Series

Colors and textures represented here are an approximate of true color under the limits of electronic viewing.

Contact ENOCHS for actual color samples. Colors subject to change without notification.

ENOCHS Treatment Tables

Horizon 2250 Horizon 49 & 49D Horizon 48



ENOCHS HORIZON

compare.









2250 High-Low Table

49D Extra-Wide Table with Drawers

49 Extra-Wide Table

49 Economy Table

Warranty:	4-YEAR			
Type:	Adjustable Height Fixed Height			
Base construction:	Fully welded all steel base	H-brace, I	aminate	All steel
Base Finish:	Powder coated	Lamir	nate	Powder coated
Top Design:		One-piece uph	nolstered top	
Height:	23" to 39"	32	n	33"
Width:	32"			24"
Length:	72"			
Drawers:	Not available on this model Two, 4" deep drawers with laminate fronts Not available or			ble on this model
Drawer Dimensions:	N/A 21.5" W x 15" D x 4" H N/A			N/A
Paper Roll Holder:	Yes			
Patient Load Capacity: (evenly distributed)	500 lbs. 300 lbs.			250 lbs.
Table Weight:	350 lbs.	139 lbs.	119 lbs.	78 lbs.
Shipping:	Ships complete no assembly required Assembly required			Assembly required
Top Color:	8 color choices			
Base Color:	Cream White or Pearl Crey Beige or Grey		Cream White or Pearl Gray	
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1			

ENOCHS HORIZON POWER 2250





New Features:

- → 500 LBS.

 PATIENT LOAD
- → High Performance Actuator System
 - ➤ Power Height
 - ➤ Smooth & Quiet
- → Foot Control
 - > Standard Feature
 - ➤ Safety Extra Low Voltage (SELV)
 - ➤ Utilizes Natural Left to Right Foot Movement

Horizon 2250 Table Also Features:

- ➤ 4-YEAR WARRANTY
- ➤ Use of Circuit Breaker for easy reset NO FUSES
- ➤ All Steel Base Construction NO PLASTIC
- One-Piece Upholstered Top, no cross seams, embossing or stitching -FOR EASY CLEANING
- ➤ 32" wide Top

ENOCHS HORIZON POWER 2250

High-Low Treatment Table

PATIENT LOAD: 500 LBS.

evenly distributed on the table's main couch





The height adjustment range of 16" allows physician or staff member to position patient at a convenient working

height. The Horizon 2250 is an excellent table for exam, procedural, treatment rooms as well as physical therapy and casting areas. The Easy Touch SELV Foot Control allows the operator to use the natural side-to-side movement of the foot to adjust the table's height from 23" to 39". Use of a foot control reduces incidents of cross contamination and enhances infection control procedures. Control pedal is covered with easy to clean non-slip tread.

TABLE FEATURES:

- 4-YEAR Warranty
- Fully welded All Steel Base Construction NO PLASTIC
- Full Power adjustment of height
- Easy to clean One-Piece Upholstered Top with 1.5" high density foam
- 32" Extra Wide Top for patient comfort
- SELV Foot Control

The Horizon 2250 High-Low Treatment Table is an ideal choice for any medical facility requiring low access for their patients. The table's adjustable height and the 32" extra wide upholstered top are two distinctive features designed for patient comfort as well as physician and staff convenience. Table has a 500 pound patient load and offers easy low access by disabled or geriatric persons in wheelchairs or with walking aids, obese patients, and those requiring physical therapy. The low 23" access height helps reduce back stress or injury to staff assisting patients on and off table.



Warranty:	4-YEAR	
Base Construction:	Fully welded all steel base	
Power System:	Electro-Mechanical provides smooth and quiet movement	
Power Adjustments:	Height: 23" to 39"	
Foot Control:	"Hands-free", standard on every table, SELV	
Hand Control:	Not available on this model	
Couch Top:	One-piece upholstered top	
Perineal Recess:	Not available on this model	
Pelvic Tilt:	Not available on this model	
Leg Extension:	Not available on this model	
Irrigation Basin:	Not available on this model	
Electrical Outlets: with Circuit Breaker	Not available on this model	
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	Not available on this model	
Stirrups:	Not available on this model	
Storage:	Not available on this model	
	Total Storage Space: N/A	
Paper Roll Holder:	Accommodates 18" or 21" paper widths	
Accessories:	Not available for this model	
Mains Power Entry: with Circuit Breaker	IEC connection, internal circuit breaker with reset	
	Table: 1 minute ON / 9 minutes OFF	
Duty Cycle:	Receptacles: N/A	
	Total receptacle load: N/A	
Ratings:	Table: 115V 60 Hz 10 amps Receptacles: N/A	
Katings.	Total receptacle load: N/A	
Table Dimensions:	Width: 32", Height: 23" - 39", Length: 72"	
Table Weight:	350 lbs.	
Patient Load:	500 lbs. evenly distributed on main couch	
Certifications:	Classified by LIL to LIL60601-1 JEC/EN 60601-1	











BASE COLORS



Colors and textures represented here are an approximate of true color under the limits of electronic viewing.

Contact ENOCHS for actual color samples. Colors subject to change without notification

ENOCHS, Inc. 800-428-2305 www.enochsmed.com enochs@enochsmed.com ©2020 48





BASE COLORS



Colors and textures represented here are an approximate of true color under the limits of electronic viewing. Contact ENOCHS for actual color samples. Colors subject to change without notification






Windsor Wine




BASE COLORS

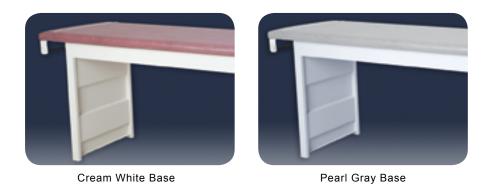

Colors and textures represented here are an approximate of true color under the limits of electronic viewing. Contact ENOCHS for actual color samples. Colors subject to change without notification

HORIZON TREATMENT Tables

Build Order Number

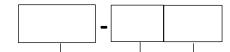

CODE	TABLE MODEL
2250	Power high-low treatment table
49D	Extra wide treatment table with drawers
49	Extra wide treatment table
48	Economy treatment table

CODE	BASE COLOR for 2250 & 48
PG	Pearl Grey
CW	Cream White

CODE BASE COLOR for 49 Series

GR Grey BE Beige

CODE	TOP COLOR
ВТ	British Tan
СВ	Cambridge Blue
EM	Emerald
FN	Fawn
II	India Ink
NN	Newberry Navy
NI	Ninja
OJ	Ocean Jade
ST	Simply Taupe
SM	Smoke
SG	Sussex Green
WW	Windsor Wine

ENOCHS Fixed-Height Tables

> Regency 456 Regency 450



ENOCHS REGENCY

compare.

	REGENCY 456	REGENCY 450	
Warranty:	4-YEAR	2-YEAR	
Base Construction:	Fully welded	all steel base	
Back Lift Mechanism:	Horizontal to 80° by spring	g assisted counter balance	
Couch Top:	100% Seamless vacuum formed s	style, all steel frame construction	
Perineal Recess:	Moulded, no seams, 10"	W x 2" D at center point	
Pelvic Tilt:	7°,	11°	
Leg Extension:	Extends out 20" fron	n end of seat section	
Irrigation Basin:	4" Dee	p basin	
Electrical Outlets:	Hospital grade, UL Listed, dual electrical outlets	Not available on this model	
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	High efficiency warmer utilizes top footend drawer	Not available on this model	
Stirrups:	Multi-positional, lateral	l and length adjustment	
Storage:	3 Side drawers: 21.5" W x 15" D x 4" H 2 Footend drawers: 15.5" W x 15" D x 4" H Total Storage Space: 3.32 cubic feet		
Foot Step:	Extra long surface engineered for 600 lb. patient load		
Paper Roll Holder:	Accommodates 18" or 21" paper widths		
Accessories:		knee crutches, padded arm board, drawer locks, stainless steel basin, restraint straps, IV pole, urology basin with screen and drain, headend accessory rails, footend accessory rails, caster base	
Mains Power Entry:	IEC connection, internal circuit breaker with reset	Not applicable	
Duty Cycle:	Receptacles: Continuous operation Instrument Warming Drawer: Continuous operation	Receptacles: Not applicable Instrument Warming Drawer: Not applicable	
Ratings:	Table: 115V 60 Hz 10 amps Receptacles: 115V 60 Hz Total receptacle load: 9 amps	Table: Not applicable Receptacles: Not applicable Total receptacle load: Not applicable	
Table Dimensions:	Width: 27", Height: 34", Length 52" (extends to 72")		
Table Weight:	268 lbs.		
Patient Load:	600 lbs. evenly distributed on main couch		
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1		



Regency 456 Table Also Features:

4-YEAR WARRANTY

- Adjustable Back with Spring Assisted Mechanism for easy operation NO GAS CYLINDER
- ➤ All Steel Base Construction NO PLASTIC
- Use of Circuit Breakers on Outlets and Warmer Drawer for easy reset NO FUSES
- ➤ 100% Seamless Vacuum Formed Top with moulded Perineal Recess
- Pelvic Tilt 7° and 11°
- Multi-Position Stirrups with length and lateral adjustment
- Pull out padded Leg Extension
- 4" deep Irrigation Basin moves independent of Leg Extension
- Five spacious Drawers with a total of 3.73 cubic feet of Storage

Pictured above: Beige Fronts, Cream White Base, Desert Sand Top

New Features:

- → 600 LBS.

 PATIENT LOAD
- → Field Reversible Side Drawers
- → Electrical Outlets
 - ➤ 9 Amp Total Load
 - ➤ Dual, Hospital Grade
- → Instrument Warming Drawer
 - Efficient Warming Technology
 - On/Off Switch with Integrated Circuit Breaker
- → Extra Long Foot Step
 - Designed for Heavier
 Patient Load

ENOCHS REGENCY 456

OB/GYN Examination Table

PATIENT LOAD: 600 LBS.

evenly distributed on the table's main couch



- 4-YEAR Warranty
- · Dual Electrical Outlets, Hospital Grade, 9 amps total load
- Instrument Warming Drawer efficiently warms drawer contents to approximate body temperature
- ON/OFF Instrument Warmer Switch with integrated Circuit Breaker - NO FUSES
- · All Steel Base and Top Construction - NO PLASTIC
- Easy to clean 100% Seamless Vacuum Formed Top with extra padding on edges for patient comfort
- . Pelvic Tilt 7" & 11"
- · Irrigation Basin designed to move independent of Leg Extension
- · Easy to adjust Stirrups with length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers



The extra long pull-out foot step features a steel stabilizing leg designed to accommodate a 600 pound patient load capacity. The step easily extends out from the table and locks into position providing the patient with a stable surface when getting on and off the table. The step's all steel construction offers extraordinary strength and durability. The 17.5' long surface is covered in easy to clean slipresistant tread.



- NO PLASTIC PARTS

New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron® on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



Located at footend of table, the high efficiency instrument warmer is designed to warm speculums and instruments to approximate body temperature. The ON/OFF switch for the warmer is located on the front side of table for instant confirmation of unit's status. An internal circuit breaker in the switch will trip if circuit is overloaded, interrupting power to the warming unit but maintaining power to the electrical receptacles.

Warranty:	4-YEAR	
Base Construction:	Fully welded all steel base	
Back Lift Mechanism:	Horizontal to 80° by spring assisted counter balance	
Couch Top:	100% Seamless vacuum formed style All steel frame construction	
Perineal Recess:	Moulded, no seams, 10" W x 2" D at center point	
Pelvic Tilt:	7°, 11°	
Leg Extension:	Extends out 20" from end of seat section	
Irrigation Basin:	4" Deep basin	
Electrical Outlets:	dual electrical outlets, hospital grade, UL Listed	
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	High efficiency warmer utilizes top footend drawer	
Stirrups:	Multi-positional, lateral and length adjustment	
Storage:	3 Side Drawers: 21.5" W x 15" D x 4" H 2 Footend Drawers: 15.5" W x 15" D x 4" H Total Storage Space: 3.32 cubic feet	
Foot Step:	Extra long surface engineered for 600 lb. patient load	
Paper Roll Holder:	Accommodates 18" or 21" paper widths	
Accessories:	knee crutches, arm board, drawer locks, stainless steel basin, restraint straps, IV pole, urology basin with drain and screen, headend accessory rails, footend accessory rails, caster base	
Mains Power Entry: with Circuit Breaker	·	
Duty Cycle:	Receptacles: Continuous Instrument Warming Drawer: Continuous	
Ratings:	Table: 115V 60Hz 10 amps Receptacles: 115V 60Hz Total receptacle load: 9 amps	
Table Dimensions:		
Table Weight:	268 lbs.	
Patient Load:	600 lbs. evenly distributed on main couch	
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1	









New Features:

- → 600 LBS.

 PATIENT LOAD
- → Field Reversible Side Drawers
- → Extra Long Foot Step
 - Designed for Heavier Patient Load
 - Easy accessibility on and off table
 - Ribbed Rubber Tread for easy cleaning

Regency 450 Table Also Features:

2-YEAR WARRANTY

- Adjustable Back with Spring Assisted Mechanism for easy operation NO GAS CYLINDER
- All Steel Base Construction NO PLASTIC
- ▶ 100% Seamless Vacuum Formed Top with moulded Perineal Recess
- Pelvic Tilt 7° and 11°
- Multi-Position Stirrups with length and lateral adjustment
- Pull out padded Leg Extension
- 4" deep Irrigation Basin moves independent of Leg Extension
- Five spacious Drawers with a total of 3.73 cubic feet of Storage

ENOCHS REGENCY 450

OB/GYN Examination Table

PATIENT LOAD: 600 LBS.

evenly distributed on the table's main couch

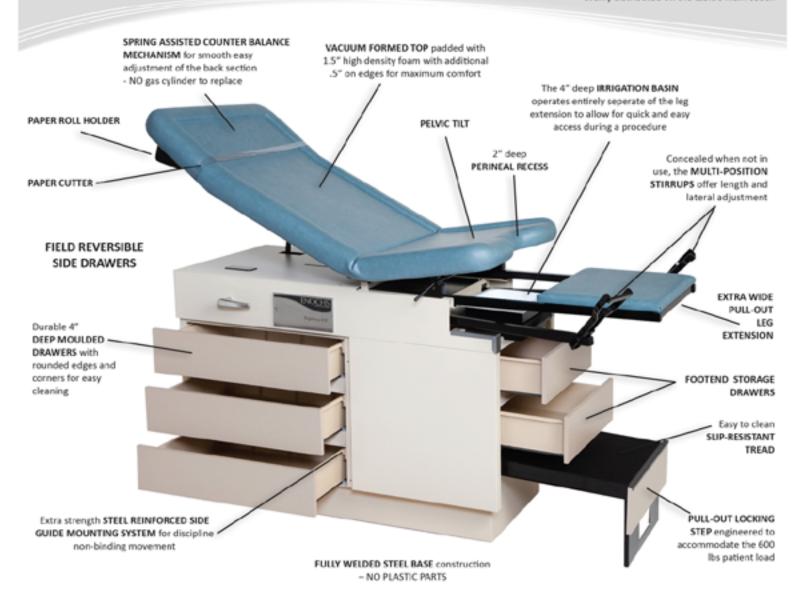


TABLE FEATURES:

- 2-YEAR Warranty
- Fully Welded all Steel Base and Top Construction - NO PLASTIC
- Spring Assisted counter balance allows smooth adjustment of couch back
- Easy to clean 100% Seamless Vacuum Formed Top
- . Pelvic Tilt with a 7" and 11" adjustment
- 4" deep Irrigation Basin engineered to move independent of Leg Extension
- Easy to adjust Stirrups with length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers



The working end of the table features the independent operation between the stirrups, leg extension and irrigation basin. This efficient design gives practitioner fast and easy transition from one to another. The padded leg extension pulls smoothly out from the table to support patient's legs and can be pushed under the seat section, allowing full access to the irrigation basin and perineal recess.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron® on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



The spring assisted counter balance lift mechanism (NO gas cylinder) allows for quick and convenient adjustment of the couch back to an infinite number of positions from horizontal to a full 80° chair position. The Vacuum formed top, designed for strength and durability has an all steel back and frame (NO wooden or plastic parts). The durable Naughaform's covered top is easy to maintain and provides a long, trouble - free service life.

Warranty:	2-YEAR	
Base Construction:	Fully welded all steel base	
Back Lift Mechanism:	Horizontal to 80° by spring assisted counter balance	
Couch Top:	100% Seamless vacuum formed style All steel frame construction	
Perineal Recess:	Moulded, no seams, 10" W x 2" D at center point	
Pelvic Tilt:	7°, 11°	
Leg Extension:	Extends out 20" from end of seat section	
Irrigation Basin:	4" Deep basin	
Electrical Outlets:	Not available on this model	
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	Not available on this model	
Stirrups:	Multi-positional, lateral and length adjustment	
Storage:	3 Side Drawers: 21.5" W x 15" D x 4" H 2 Footend Drawers: 15.5" W x 15" D x 4" H Total Storage Space: 3.32 cubic feet	
Foot Step:	Extra long surface engineered for 600 lb. patient load	
Paper Roll Holder:	Accommodates 18" or 21" paper widths	
Accessories:	knee crutches, arm board, drawer locks, stainless steel basin, restraint straps, IV pole, urology basin with drain and screen, headend accessory rails, footend accessory rails, caster base	
Mains Power Entry: with Circuit Breaker	Not available on this model	
Duty Cycle:	Receptacles: N/A Instrument Warming Drawer: N/A	
Ratings:	Table: N/A Receptacles: N/A Total receptacle load: N/A	
Table Dimensions:	Width: 27", Height: 34", Length 52" (extends to 72")	
Table Weight:	268 lbs.	
Patient Load:	600 lbs. evenly distributed on main couch	
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1	








REGENCY & EXCEL Accessories

Model #6109 Universal Clamp



Model #6110 Tri-Clamp

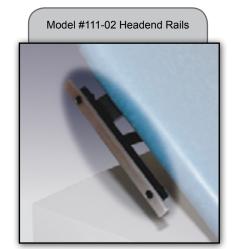

Model #111-03 Footend Rails



Factory Installed Only




For Regency Table Only Requires #111-02 Headend Rails and #6109 or #6110 Clamp

Model #111-10 Urology Pan

For Regency Table Only



Model #111-11 Knee Crutches

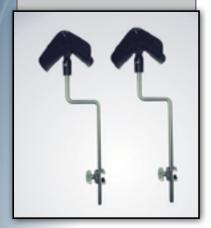

Requires #111-03 Footend Rails

Model #111-12 Restraint Straps (pair)



60





Cambridge Blue



Smoke



Cambridge Blue (above)



Windsor Wine



English Brown



Desert Sand



Smoke (above)



Ocean Jade



Highland Sage



India Ink



Newberry Navy



Simply Taupe

BASE/DRAWER FRONT COLORS



Cream White Base Beige Drawer Fronts



Pearl Gray Base Gray Drawer Fronts




Cream White Base Natural Oak Drawer Fronts

Colors and textures represented here are an approximate of true color under the limits of electronic viewing.

Contact ENOCHS for actual color samples. Colors subject to change without notification

REGENCY OB/GYN Examination Table

Build Order Number

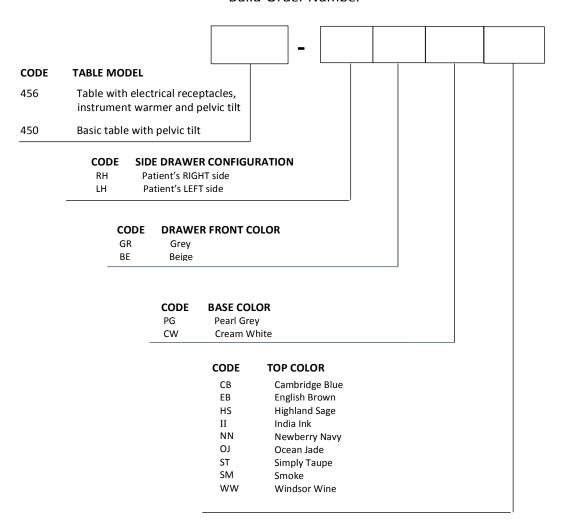

ENOCHS Fixed-Height Tables

> Excel 356 Excel 355 Excel 350



ENOCHS EXCEL

compare.

	EXCEL 356	EXCEL 355	EXCEL 350	
Warranty:	4-YEAR	4-YEAR	2-YEAR	
Base construction:		Fully welded all steel base		
Back Lift Mechanism:	Но	rizontal to 80° by spring assisted counter balan	ce	
Couch Top:		One-piece upholstered top, steel reinforced		
Perineal Recess:		Not available on this model		
Pelvic Tilt:	11°	Not available	on this model	
Leg Extension:		Extends out 18" from end of seat section		
Irrigation Basin:		4" Deep basin		
Electrical Outlets:	Hospital grade, UL Liste	d, dual electrical outlets	Not available on this model	
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	High efficiency warmer utilizes top footend drawer Not available on this model		on this model	
Stirrups:		Multi-positional, lateral and length adjustment		
Storage:		3 Side drawers: 21.5" W x 15" D x 4" H 2 Footend drawers: 15.5" W x 15" D x 4" H Total Storage Space: 3.32 cubic feet		
Foot Step:	Ext	Extra long surface engineered for 600 lb. patient load		
Paper Roll Holder:		Accommodates 18" or 21" paper widths		
Accessories:	knee crutches, drawer locks, stainless steel b	knee crutches, drawer locks, stainless steel basin, restraint straps, IV pole, urology basin with drain and screen, footend accessory rails, caster base		
Mains Power Entry:	IEC connection, internal circuit breaker with reset	IEC connection, internal circuit breaker with reset	Not applicable	
Duty Cycle:	Receptacles: Continuous operation Instrument Warming Drawer: Continuous operation	Receptacles: Continuous operation Instrument Warming Drawer: Not applicable	Receptacles: Not applicable Instrument Warming Drawer: Not applicable	
Ratings:	Table: 115V 60 Hz 10 ampsTable: Not applicableReceptacles: 115V 60 HzReceptacles: Not applicableTotal receptacle load: 9 ampsTotal load: Not applicable		Receptacles: Not applicable	
Table Dimensions:	Width: 27", Height: 34", Length 55" (extends to 73")			
Table Weight:		268 lbs.		
Patient Load:		600 lbs. evenly distributed on main couch		
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1			



Excel 356 Table Also Features:

- ▶ 4-YEAR WARRANTY
- Adjustable Back with Spring Assisted Mechanism for easy operation NO GAS CYLINDER
- Use of Circuit Breakers on Outlets and Warmer Drawer for easy reset NO FUSES
- All Steel Base Construction NO PLASTIC
- Easy to clean One-Piece Upholstered Top designed without cross-seams, decorative embossing or stitching
- ➤ Pelvic Tilt 11°
- Multi-Position Stirrups with length and lateral adjustment
- Pull out padded Leg Extension
- 4" deep Irrigation Basin moves independent of Leg Extension
- Five spacious Drawers with a total of 3.73 cubic feet of Storage

Pictured above: Beige Fronts, Cream White Base, Ninia Top

New Features:

- → 600 LBS.

 PATIENT LOAD
- → Field Reversible Side Drawers
- → Electrical Outlets
 - > 9 Amp Total Load
 - > Dual, Hospital Grade
- → Instrument Warming Drawers
 - Efficient Warming Technology
 - On/Off Switch with Integrated Circuit Breaker
- → Extra Long Foot Step
 - Designed for Heavier Patient Load

ENOCHS EXCEL 356

General Examination Table

PATIENT LOAD: 600 LBS.

evenly distributed on the table's main couch



- · 4-YEAR Warranty
- · All welded Steel Base treated for resistance to rust
- . Dual Electrical Outlets, Hospital Grade, 9 amps total load
- . Instrument Warming Drawer efficiently warms drawer contents to approximate body temperature
- · ON/OFF Warmer Switch with integrated Circuit Breaker - NO FUSES
- · Easy to clean One-Piece Upholstered Top designed without cross-seams, decorative stitching or embossing
- 11° Pelvic Tift with release levers on both sides of table
- · Easy to adjust Stirrups with length and lateral adjustment
- · 3.32 cubic feet of Storage in the five spacious Drawers



The extra long pull-out foot step features a steel stabilizing leg designed to accommodate a 600 pound patient load capacity. The step easily extends out from the table and locks into position providing the patient with a stable surface when getting on and off the table. The step's all steel construction offers extraordinary strength and durability. The 17.5" long surface is covered in easy to clean slipresistant tread.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two foctend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron* on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



The spring assisted counter balance lift mechanism (NO gas cylinder) allows for quick and convenient adjustment of the couch back to an infinite number of positions from horizontal to a full 80° chair position. The table's couch top is designed for strength and durability with a steel reinforced frame and covered in durable Naughahyde* material that provides a long, trouble - free service

Specifications for the ENOCHS EXCEL 356 General Examination Table

Warranty:	4-YEAR	
Base Construction:	Fully welded all steel base	
Back Lift Mechanism:	Horizontal to 80° by spring assisted counter balance	
Couch Top:	One-piece upholstered top, steel reinforced	
Perineal Recess:	Not available on this model	
Pelvic Tilt:	11°	
Leg Extension:	Extends out 18" from end of seat section	
Irrigation Basin:	4" Deep basin	
Electrical Outlets:	Dual electrical outlets, hospital grade, UL Listed	
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	High efficiency warmer utilizes top footend drawer	
Stirrups:	Multi-positional, lateral and length adjustment	
Storage:	3 Side Drawers: 21.5" W x 15" D x 4" H 2 Footend Drawers: 15.5" W x 15" D x 4" H Total Storage Space: 3.32 cubic feet	
Foot Step:	Extra long surface engineered for 600 lb. patient load	
Paper Roll Holder:	Accommodates 18" or 21" paper widths	
Accessories:	knee crutches, drawer locks, stainless steel basin, restraint straps, IV pole, urology basin with drain and screen, footend accessory rails, caster base	
Mains Power Entry: with Circuit Breaker	Standard feature	
Duty Cycle:	Receptacles: Continuous Instrument Warming Drawer: Continuous	
Ratings:	Table: 115V 60Hz 10 amps Receptacles: 115V 60Hz Total receptacle load: 9 amps	
Table Dimensions:	Width: 27", Height: 34", Length 55" (extends to 73")	
Table Weight:	268 lbs.	
Patient Load:	600 lbs. evenly distributed on main couch	
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1	







ENOCHS EXAMINING ROOM FURNITURE



New Features:

- → 600 LBS.

 PATIENT LOAD
- → Field Reversible Side Drawers
- → Electrical Outlets
 - ➤ 9 Amp Total Load
 - ➤ Dual, Hospital Grade
- → Extra Long Foot Step
 - Designed for Heavier
 Patient Load
 - ➤ Easy Access on and off table
 - ➤ Ribbed Rubber Tread for easy cleaning

Excel 355 Table Also Features:

- 4-YEAR WARRANTY
- Adjustable Back with Spring Assisted Mechanism for easy operation NO GAS CYLINDER
- Use of Circuit Breakers for easy reset NO FUSES
- All Steel Base Construction NO PLASTIC
- Easy to clean One-Piece Upholstered Top designed without cross-seams, decorative embossing or stitching
- Multi-Position Stirrups with length and lateral adjustment
- ➤ Pull out padded Leg Extension
- 4" deep Irrigation Basin moves independent of Leg Extension
- Five spacious drawers with a total of 3.73 cubic feet of Storage

Pictured above: Grey Fronts, Pearl Grey Base, Smoke Top

ENOCHS EXCEL 355

General Examination Table

PATIENT LOAD: 600 LBS.

evenly distributed on the table's main couch

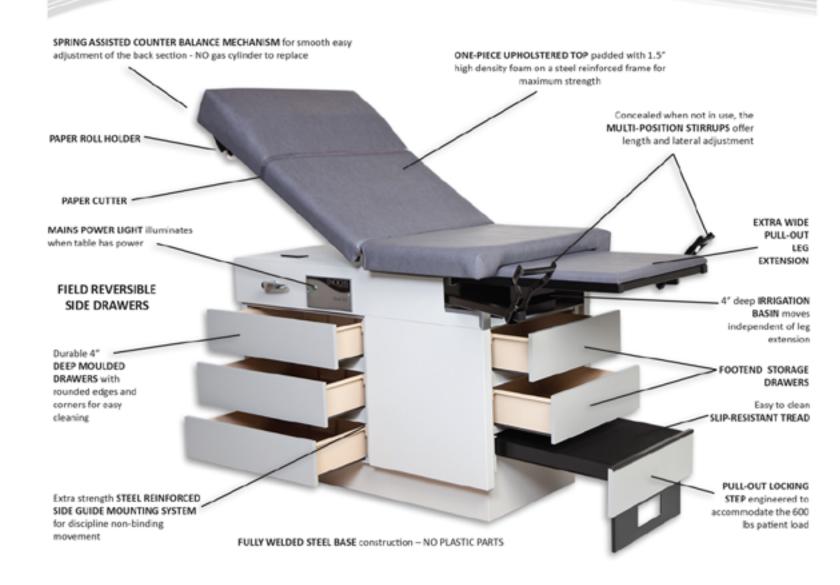


TABLE FEATURES:

- 4-YEAR Warranty
- All Steel Base Construction NO PLASTIC
- Dual Electrical Outlets, Hospital Grade, 9 amps total load
- Easy to clean One-Piece Upholstered Top, designed without cross-seams, decorative stitching or embossing
- Top constructed with a Steel Reinforced
 Frame for maximum strength
- Irrigation Basin moves independent of Leg Extension for fast transition
- Multi-position Stirrups with easy length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers



The extra long pull-out foot step features a steel stabilizing leg designed to accommodate a 600 pound patient load cepacity. The step easily extends out from the table and locks into position providing the patient with a stable surface when getting on and off the table. The step's all steel construction offers extraordinary strength and durability. The 17.5" long surface is covered in easy to clean slip-resistant tread.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron* on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



The table's electrical receptacles feature a new tool required access door designed to restrict unauthorized use of the outlets. Door easily unlocks so accessories can be plugged into the UL Listed, hospital grade receptacles. Once in place, the door closes over the electrical cords, securing plugs in place with the locking mechanism. Receptacles have a 9 amps total maximum load.

Specifications for the ENOCHS EXCEL 355 General Examination Table

Warranty:	4-YEAR	
Base Construction:	Fully welded all steel base	
Back Lift Mechanism:	Horizontal to 80° by spring assisted counter balance	
Couch Top:	One-piece upholstered top, steel reinforced	
Perineal Recess:	Not available on this model	
Pelvic Tilt:	Not available on this model	
Leg Extension:	Extends out 18" from end of seat section	
Irrigation Basin:	4" Deep basin	
Electrical Outlets:	Dual electrical outlets, hospital grade	
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	Not available on this model	
Stirrups:	Multi-positional, lateral and length adjustment	
Storage:	3 Side Drawers: 21.5" W x 15" D x 4" H 2 Footend Drawers: 15.5" W x 15" D x 4" H Total Storage Space: 3.32 cubic feet	
Foot Step:	Extra long surface engineered for 600 lb. patient load	
Paper Roll Holder:	Accommodates 18" or 21" paper widths	
Accessories:	knee crutches, drawer locks, stainless steel basin, restraint straps, IV pole, urology basin with drain and screen, footend accessory rails, caster base	
Mains Power Entry: with Circuit Breaker	Standard feature	
Duty Cycle:	Receptacles: Continuous Instrument Warming Drawer: N/A	
Ratings:	Table: 115V 60Hz 10 amps Receptacles: 115V 60Hz Total receptacle load: 9 amps	
Table Dimensions:	Width: 27", Height: 34", Length 55" (extends to 73")	
Table Weight:	268 lbs.	
Patient Load:	600 lbs. evenly distributed on main couch	
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1	









New Features:

- → 600 LBS.

 PATIENT LOAD
- → Field Reversible Side Drawers
- → Extra Long Foot Step
 - Designed for Heavier
 Patient Load
 - ➤ Easy Access on and off table
 - Ribbed Rubber Tread for easy cleaning

Excel 350 Table Also Features:

- 2-YEAR WARRANTY
- Adjustable Back with Spring Asisted Mechanism for easy operation NO GAS CYLINDER
- All Steel Base Construction NO PLASTIC
- Easy to clean One-Piece Upholstered Top designed without cross-seams, decorative embossing or stitching
- Multi-Position Stirrups with length and lateral adjustment
- Pull out padded Leg Extension
- ➤ 4" deep Irrigation Basin moves independent of Leg Extension
- Five spacious Drawers with a total of 3.73 cubic feet of Storage

ENOCHS EXCEL 350

General Examination Table

PATIENT LOAD: 600 LBS.

evenly distributed on the table's main couch

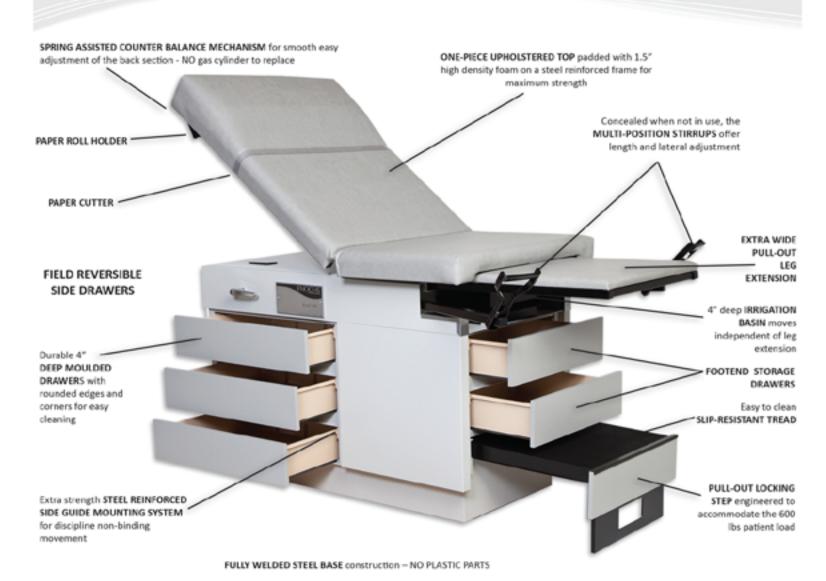


TABLE FEATURES:

- · 2-YEAR Warranty
- Pull-Out Locking Step engineered to accommodate the 600 lbs patient load
- All Steel Base Construction
 NO PLASTIC
- Spring Assisted counter balance allows smooth adjustment of couch back
- Easy to clean One-Piece
 Upholstered Top designed without cross-seams, decorative stitching or embossing
- Multi-position Stirrups with easy length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers



The extra long pull-out foot step features a steel stabilizing leg designed to accommodate a 600 pound patient load capacity. The step easily extends out from the table and locks into position providing the patient with a stable surface when getting on and off the table. The step's all steel construction offers extraordinary strength and durability. The 17.5" long surface is covered in easy to clean slip-resistant tread.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron® on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



The working end of the table has an exclusive ENOCHS design of a pullout padded leg extension that moves independent of the 4" deep irrigation basin. This separate basin design allows physician to easily transition from leg extension to basin and provides a fast easy clean up by the staff after the procedure. The Excel 350 also features multi-position stirrups that lock into position with the weight of the patient's foot.

Specifications for the ENOCHS EXCEL 350 General Examination Table

Warranty:	2-YEAR	
Base Construction:	Fully welded all steel base	
Back Lift Mechanism:	Horizontal to 80° by spring assisted counter balance	
Couch Top:	100% Seamless vacuuOne-piece upholstered top, steel reinforced	
Perineal Recess:	Not available on this model	
Pelvic Tilt:	Not available on this model	
Leg Extension:	Extends out 18" from end of seat section	
Irrigation Basin:	4" Deep basin	
Electrical Outlets:	Not available on this model	
Instrument Warmer: with ON/OFF Switch, Circuit Breaker & Reset	Not available on this model	
Stirrups:	Multi-positional, lateral and length adjustment	
Storage:	3 Side Drawers: 21.5" W x 15" D x 4" H 2 Footend Drawers: 15.5" W x 15" D x 4" H Total Storage Space: 3.32 cubic feet	
Foot Step:		
Paper Roll Holder:	Accommodates 18" or 21" paper widths	
Accessories:	knee crutches, drawer locks, stainless steel basin, restraint straps, IV pole, urology basin with drain and screen, footend accessory rails, caster base	
Mains Power Entry: with Circuit Breaker	Not available on this model	
Duty Cycle:	Receptacles: N/A Instrument Warming Drawer: N/A	
Ratings:	Table: N/A Receptacles: N/A Total receptacle load: N/A	
Table Dimensions:	Width: 27", Height: 34", Length 55" (extends to 73")	
Table Weight:	268 lbs.	
Patient Load:	600 lbs. evenly distributed on main couch	
Certifications:	Classified by UL to UL60601-1, IEC/EN 60601-1, CAN/CSA-C22.2 No. 601.1	







REGENCY & EXCEL Accessories




Model #6109 Universal Clamp



Model #6110 Tri-Clamp



Model #111-02 Headend Rails

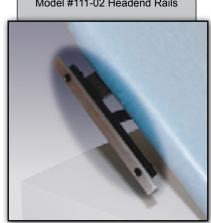

Model #111-03 Footend Rails



Factory Installed Only

Model #111-05 Arm Board




For Regency Table Only Requires #111-02 Headend Rails and #6109 or #6110 Clamp

For Regency Table Only



Model #111-11 Knee Crutches

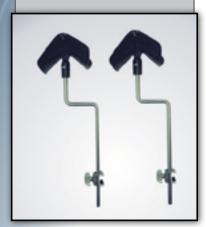

Requires #111-03 Footend Rails

Model #111-12 Restraint Straps (pair)



Model #111-50 Caster Base





EXCEL Color Selections

UPHOLSTERY COLORS





Smoke

Fawn



Cambridge Blue



Windsor Wine



Emerald



Fawn (above)



Smoke (above)



Ocean Jade



British Tan



Ninja

Sussex Green



India Ink



Newberry Navy



Simply Taupe

BASE/DRAWER FRONT COLORS



Cream White Base Beige Drawer Fronts



Pearl Gray Base Gray Drawer Fronts




Cream White Base Natural Oak Drawer Fronts

Colors and textures represented here are an approximate of true color under the limits of electronic viewing.

Contact ENOCHS for actual color samples. Colors subject to change without notification

EXCEL General Examination Table

Build Order Number

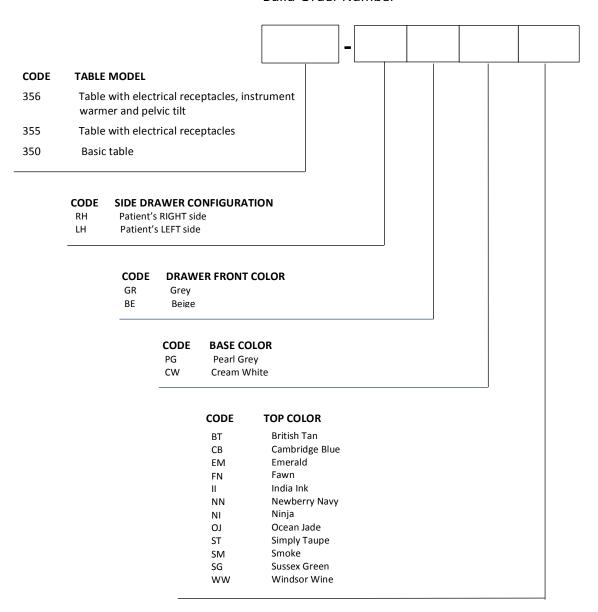

ENOCHS

ENOCHS Physician Seating

Specifications for Physician Seating









	92 OPERATOR STOOL	96 OPERATOR STOOL WITH BACK	101 PNUEMATIC STOOL	106 PNUEMATIC STOOL WITH BACK
Warranty:	4 year	4 year	4 year	4 year
Base Construction:	fully welded, tubular steel base	fully welded, tubular steel base	cast aluminum	cast aluminum
Base Finish:	durable powder coated finish	durable powder coated finish	polished aluminum	polished aluminum
Number of Legs:	5	5	5	5
Foot Ring:	yes	yes	no*	no*
Caster Style:	2" ball bearing, hard surface caster	2" ball bearing, hard surface caster	2" carpet caster	2" carpet caster
Adjustment System:	threaded stem	threaded stem	pneumatic gas cylinder	pneumatic gas cylinder
Adjustment By:	rotation of seat	rotation of seat	hand lever under seat	hand lever under seat
Height:	18.5" to 25"	18.5" to 25"	20" to 30"	20" to 30"
Back Support:	no	yes	no	yes
Back Adjustment:	N/A	3" vertical	N/A	3" vertical
Seat Dimensions:	14"W x 2.5"H	14"W x 2.5"H	14"W x 2.5"H	14"W x 2.5"H
Weight Capacity:	350 lbs.	225 lbs.	350 lbs.	225 lbs.
Stool Weight:	16 lbs.	21 lbs.	13 lbs.	18 lbs.
Shipping:	assembly required	assembly required	assembly required	assembly required
Color Selection:	Seat: Black, Cambridge Blue, Desert Sand, Earthen Brown, English Brown, Emerald, Fawn, Ninja, Ocean Jade, Smoke, Windsor Wine			
*available as an add on accessory				





PHYSICIAN SEATING Color Selections

UPHOLSTERY COLORS





India Ink



Cambridge Blue



Windsor Wine



Ninja



English Brown



Emerald



Fawn (above)



Smoke



Ocean Jade



British Tan



Sussex Green



India Ink (above)



Newberry Navy



Simply Taupe



Highland Šage

© 2020 79

Colors and textures represented here are an approximate of true color under the limits of electronic viewing. Contact ENOCHS for actual color samples. Colors subject to change without notification

ENOCHS, Inc. 800-428-2305 enochs@enochsmed.com www.enochsmed.com

PHYSICIAN SEATING

Build Order Number



CODE	Model
92	Basic 5-leg operator stool
96	Basic 5-leg operator stool with back
101	Pneumatic 5-leg operator stool
101	Pneumatic 5-leg operator stool with back

CODE	TOP COLOR
ВТ	British Tan
СВ	Cambridge Blue
EM	Emerald
EB	English Brown
FN	Fawn
HS	Highland Sage
II	India Ink
NN	Newberry Navy
NI	Ninja
OJ	Ocean Jade
ST	Simply Taupe
SM	Smoke
SG	Sussex Green
WW	Windsor Wine

ENOCHS

ENOCHS Warranty Sales Policy

The ENOCHS Guarantee of Quality

ENOCHS, Inc. guarantees ENOCHS products to be free from defects in material or workmanship under normal installation, use and service, for a period of **4 years** (2 years on selected products) from date of shipment from our factory.

All items and component parts manufactured for ENOCHS by others are warranted against defects for a period equal to warranty offered by original equipment manufacturer. This guarantee is invalid if product is damaged as a result of misuse, abuse, neglect or alteration. Transportation costs for return of product or part is responsibility of purchaser. All other expressed or implied warranties of the company are superseded.

ENOCHS Examining Room Furniture www.enochsmed.com

Mailing Address

245 Erie Street

Huntington, IN 46750

Customer Service

800-428-2305

FAX: 317-580-2944

customerservice@enochsmed.com

Technical, Service and Repair Support

800-322-6416

FAX: 260-359-2004

tech@enochsmed.com

ENOCHS WARRANTY | www.enochsmed.com/warranty

ENOCHS, Inc. guarantees ENOCHS products to be free from defects in material or workmanship under normal installation, use and service, for a period of four years from date of shipment from our factory. All items and component parts manufactured for ENOCHS by others are warranted against defects for a period equal to warranty offered by original equipment manufacturer. This guarantee is invalid if product is damaged as a result of misuse, abuse, neglect or alteration. Transportation costs for return of product or part is responsibility or purchaser. All other expressed or implied warranties of the company are superseded.

ENOCHS Sales Policy | www.enochsmed.com/policy

ENOCHS sales are made through authorized medical supply distributors. Call ENOCHS Customer Service for Distributor for your area. ENOCHS reserves the right to bid direct to state and federal agencies.

Ordering Information

Distributor Purchase Order to ENOCHS must include ENOCHS catalog number, designation of side drawers on patient's right or left, color of drawer fronts, base and top where applicable. Orders must include Consignee name, address, delivery contact with phone number and any other special shipping instructions. Orders may be e-mailed or FAXED to ENOCHS. Phone orders are not accepted.

Custom Orders

Products ordered with special custom upholstery and/or drawer fronts will incur additional charges and be subject to delay in delivery. Sample of each custom color is required at time order is placed. Such orders, once accepted and in process, are not subject to cancellation.

Special Products

Orders for furniture of special construction or finish are not subject to cancellation after acceptance of ENOCHS acknowledgment.

Cancellations

Standard orders canceled five (5) or more days after entry of order will be subject to a 25% cancellation fee.

Freight Charges

All shipments will be freight collect except when a drop shipment to a third party requested. All bill of lading drop shipments directly to enduser will be marked as third party billing, allowing delivering carrier to collect freight charges directly from dealer whose name and address appear on bill of lading. This policy will be followed unless otherwise negotiated.

Export Packaging

Export packaging is available at an additional fee.

Transportation Damages

All ENOCHS products are securely packaged for safe transit. ENOCHS holds a clear bill of lading when shipment leaves our dock. ENOCHS must be notified of the damage by Consignee or distributor within 5 days of receipt of freight. ENOCHS will execute freight claims on goods shipped on a prepaid basis. Freight claims on collect or third party billing shipments will be the full responsibility of Consignee or distributor.

Visible Damage Upon delivery, Consignee is responsible to check outside carton of each item for any visible damage before accepting an ENOCHS shipment. Carton damage might indicate merchandise damage – BEFORE ACCEPTING SHIPMENT MAKE NOTATION ON DELIVERY RECEIPT OF THE VISIBLE DAMAGE. Call ENOCHS Customer Service immediately to register damage.

Concealed Damage

Upon delivery, Consignee is responsible for inspecting outside of carton of each item for any visible damage. If carton damage indicates possible merchandise damage, note condition of carton on delivery receipt. UNPACK merchandise immediately and check for damage. If damage is discovered, it is the receiver's responsibility to follow procedure to insure proper credit and filing of claim.

- 1. Call delivering carrier for an immediate freight inspection.
- 2. Do not move merchandise and retain all packing material.
- 3. Call ENOCHS immediately to register concealed damage and for further instruction.

ENOCHS will accept damaged merchandise for estimate of repair for freight claim only if steps 1-3 have been followed.

Return Merchandise Policy and Procedure

ENOCHS will accept merchandise for credit review if procedure is correctly followed. All return merchandise is subject to a 25% restocking charge. Credit will reflect merchandise price, less restocking fee and all charges for repairs or replacements. ENOCHS reserves the right to refuse returned merchandise if procedure is not explicitly followed.

- 1. Call ENOCHS within 30 days of original ship date for return authorization number. Returns without proper authorization will be refused.
- 2. Returning party is responsible for packing merchandise for safe transit. All damage due to improper packing will be deducted from credit.
- 3. Merchandise must be received by ENOCHS within 30 days from date of return authorization.
- 4. All freight charges must be prepaid.

ENOCHS Product Design

ENOCHS reserves the right to make product improvements or change product specification without notification.

Prices

Prices are effective on date of shipment and not on date entry of order. All prices are F.O.B. Factory and are subject to change without notice.

ENOCHS

click here to

LOCATE ITEM ON WEBSITE

ENOCHS Examining Room Furniture www.enochsmed.com

Mailing Address

245 Erie Street Huntington, IN 46750

Customer Service

800-428-2305

FAX: 317-580-2944

customerservice@enochsmed.com

Technical, Service and Repair Support

800-322-6416

FAX: 260-359-2004

tech@enochsmed.com



www.facebook.com/enochsmed