
## **ENOCHS ENCORE 4200**

## **High-Low Examination Table**

PATIENT LOAD: 500 LBS.

evenly distributed on the table's main couch

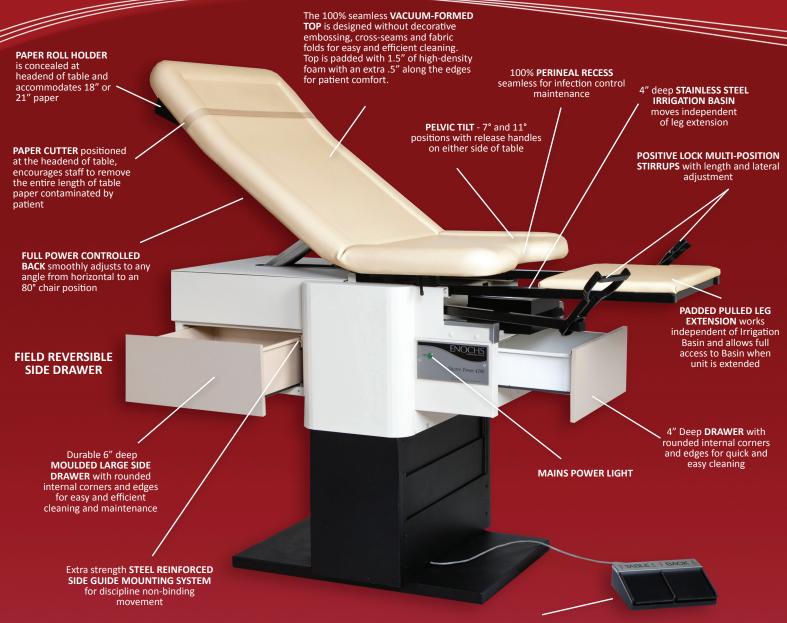

**'EASY TOUCH' FOOT CONTROL** adjusts table's height and back with reliable hands-free action. Pedals are covered with non-slip tread and utilize a natural left to right rocking motion of operator's foot.

## **TABLE FEATURES:**

- 2-YEAR Warranty
- Hand Control available as an accessory
- Fully welded All Steel Base Construction NO PLASTIC
- Full Power adjustment of back
- Easy to clean 100% Seamless Vacuum Formed Top with additional padding on edges for patient comfort
- Pelvic Tilt 7° and 11°
- Stainless Steel Irrigation Basin operates independent of Leg Extension
- Easy to adjust Stirrups with length and lateral adjustment
- 1.4 cubic feet of Storage



Securely attached to the table base, the ENOCHS vacuum-formed top is the strongest top on the market with an all steel frame and back. Top is padded for patient comfort with 1.5" of high density foam with an additional .5" of foam contouring the edges. The sleek Euro-look top features a seamless surface, moulded perineal recess and a pelvic tilt with positions of 7° and 11°. The top offers proven performance of strength, durability and long service life.



The working end of the table has the exclusive ENOCHS design of the pullout padded leg extension moving independent of the 4" deep stainless steel irrigation basin. This design reduces surface area to allow fast easy clean up after irrigation procedure. Patients with limited rotation can be placed in a GYN position utilizing the pad board instead of stirrups while physician maintains unrestricted access to the irrigation basin.