

ENOCHS ENCORE 4250

High-Low Examination Table

PATIENT LOAD: 500 LBS.
evenly distributed on the table's main couch


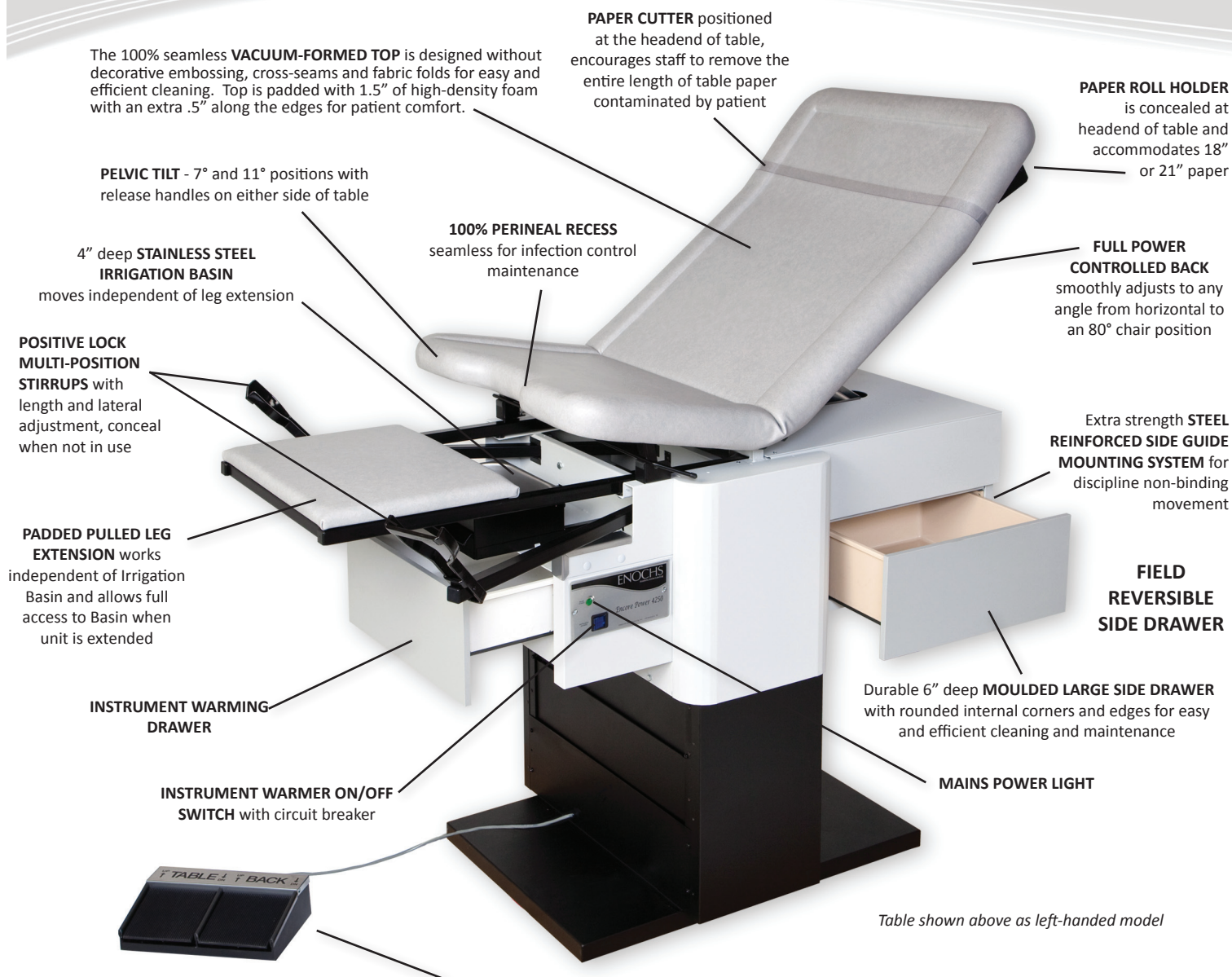


Table shown above as left-handed model

'EASY TOUCH' FOOT CONTROL adjusts table's height and back with reliable hands-free action. Pedals are covered with non-slip tread and utilize a natural left to right rocking motion of operator's foot.

TABLE FEATURES:

- 4-YEAR Warranty
- Hand Control available as an accessory
- Fully welded All Steel Base Construction - NO PLASTIC
- Full Power adjustment of back
- Easy to clean 100% Seamless Vacuum Formed Top with additional padding on edges for patient comfort
- Pelvic Tilt 7° and 11°
- Electrical Outlets with circuit breaker - NO FUSES
- Instrument Warming Drawer efficiently warms drawer contents to approximate body temperature
- Instrument Warmer ON/OFF Switch with Integrated Circuit Breaker - NO FUSES
- Stainless Steel Irrigation Basin operates independent of Leg Extension
- Easy to adjust Stirrups with length and lateral adjustment
- 1.4 cubic feet of Storage

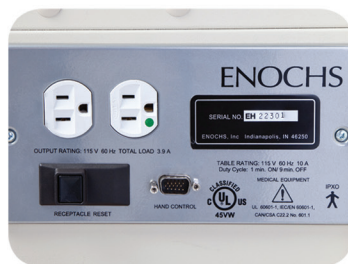

The hospital grade dual electrical outlets are conveniently located at the foot end of the table. The outlets feature a dedicated circuit breaker that will trip if outlets are overloaded. A quick reset of the circuit breaker re-establishes power to the outlets, eliminating costly down time. The Power 4250 also features a SELV (safety extra low voltage) hand control, an accessory that can be ordered at time of purchase or added at a later date.



The footend of the table features a padded pullout leg extension, 4\"/>