


ENOCHS EXCEL 350

General Examination Table

PATIENT LOAD: 600 LBS.
evenly distributed on the table's main couch

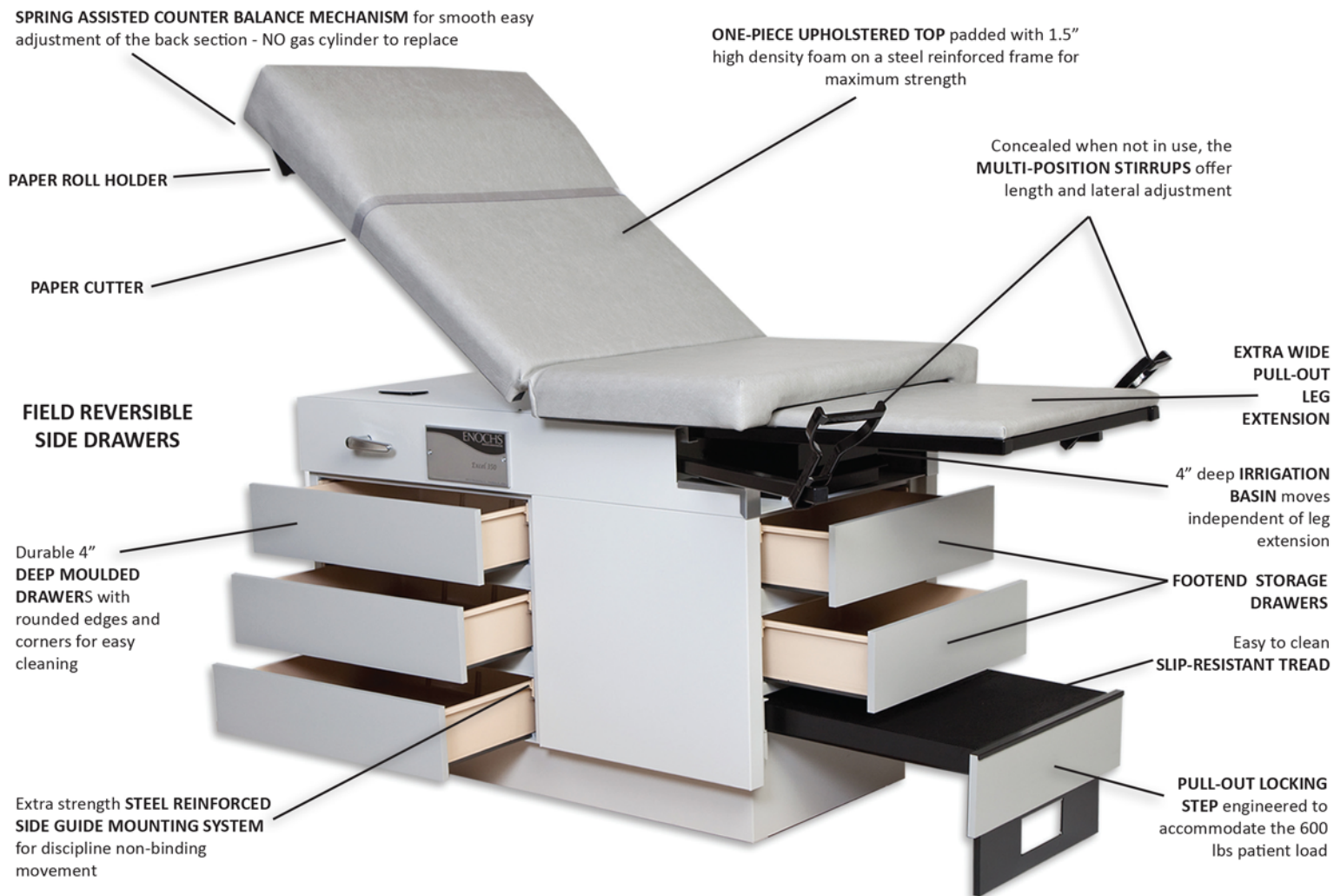

FULLY WELDED STEEL BASE construction – NO PLASTIC PARTS

TABLE FEATURES:

- 2-YEAR Warranty
- Pull-Out Locking Step engineered to accommodate the 600 lbs patient load
- All Steel Base Construction – NO PLASTIC
- Spring Assisted counter balance allows smooth adjustment of couch back
- Easy to clean One-Piece Upholstered Top designed without cross-seams, decorative stitching or embossing
- Multi-position Stirrups with easy length and lateral adjustment
- 3.32 cubic feet of Storage in the five spacious Drawers



The extra long pull-out foot step features a steel stabilizing leg designed to accommodate a 600 pound patient load capacity. The step easily extends out from the table and locks into position providing the patient with a stable surface when getting on and off the table. The step's all steel construction offers extraordinary strength and durability. The 17.5" long surface is covered in easy to clean slip-resistant tread.



New field reversible side drawers allows table to have the flexibility to be placed into any room configuration. Side drawers combined with the two footend drawers offer maximum storage capacity totaling 3.32 cubic feet. Nylatron® on the bottom drawer edges protect against wear and enhances the drawer's overall performance. An integrated stops eliminate accidental spillage of drawer contents.



The working end of the table has an exclusive ENOCHS design of a pullout padded leg extension that moves independent of the 4" deep irrigation basin. This separate basin design allows physician to easily transition from leg extension to basin and provides a fast easy clean up by the staff after the procedure. The Excel 350 also features multi-position stirrups that lock into position with the weight of the patient's foot.